


Glycols in Toothpaste Standard on ZB-WAXplus**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:** 130 °C to 220 °C @ 20 °C/min for 2 min.**Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Split 20:1 0.5 μ L @ 225°C**Detection:** Mass Spectrometer (MS) (180°C)**Analyst Note:** 2.5 mg/ml glycols standards plus 0.5 gm/ml internal standard.
added 5 m retention gap to column.**Products used in this application:**

16507

**ANALYTES:**

- 1** 1,2-Propanediol
- 2** Ethylene glycol (EG)
- 3** 1,3-Propanediol
- 4** Diethyleneglycol
- 5** Glycerin

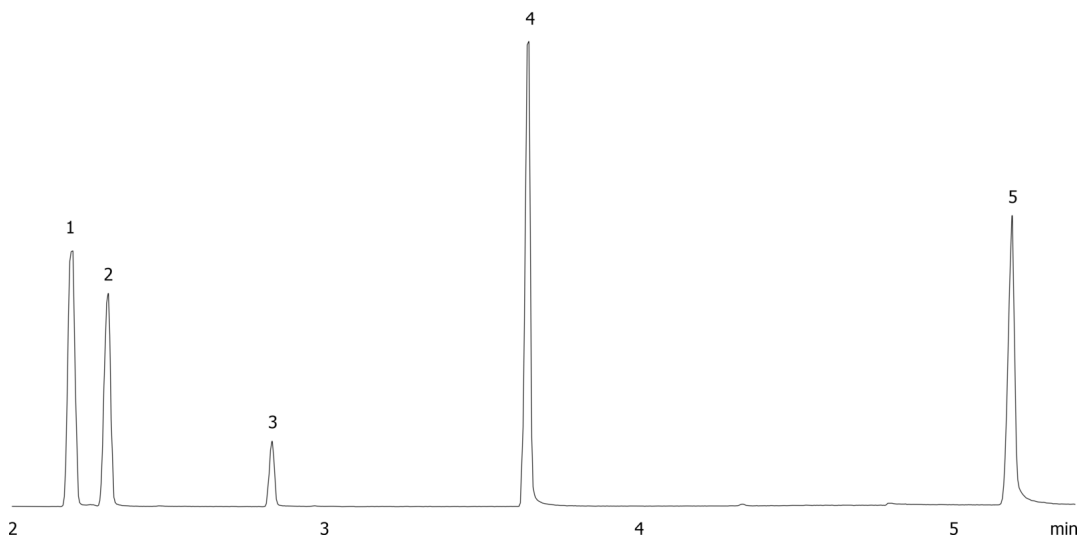