


Sterols from Olive Oil on ZB-MR2**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 230 °C to 320 °C @ 10 °C/min for 4 min**Carrier Gas:** Constant Flow Helium, 1.3 mL/min**Injection:** Splitless :1 0.5 μ L @ 300°C**Detection:** Mass Spectrometer (MS) (180°C)**Analyst Note:** Sample prepared using Saponification, Solid Phase Extraction (SPE), and Derivatization.**Products used in this application:**

16543

**ANALYTES:**

- 1 Stigmasterol
- 2 b-Sitosterol
- 3 Betulin

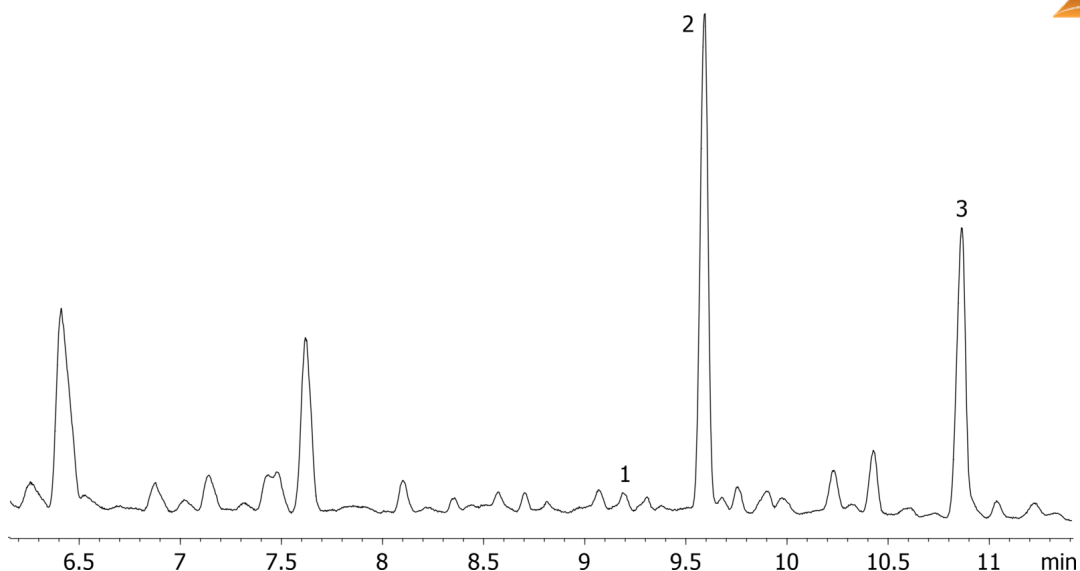