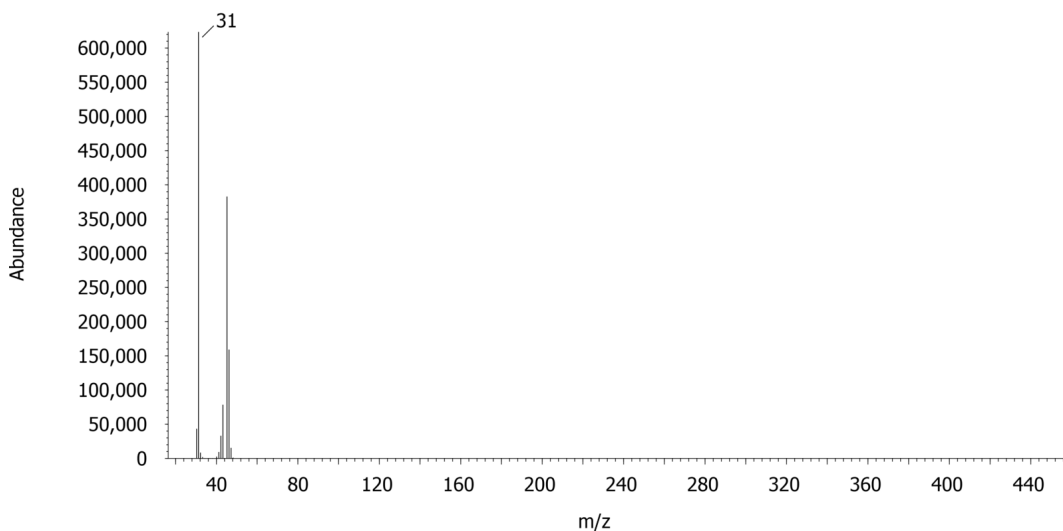


Ethanol by MSD on Zebron ZB-Bioethanol Column**Column:** Zebron™ ZB-Bioethanol, GC Cap. Column 30 m x 0.25 mm x 1.00 µm, Ea**Phase:** Proprietary Bioethanol Phase**Dimensions:** 30 meters x 0.25 mm x 1 µm**Order No:** 7HG-G020-22**Oven Profile:** 40°C for 5 min to 300°C @ 25°C/min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 100:1 0.1 µL @ 240°C**Detection:** Mass Selective (MSD) (230°C)**Analyst Note:** Ethanol ~95%, Methanol ~2%, Gasoline ~2%, Fusel Alcohols ~1%

16674

**ANALYTES:****1** Ethanol**Products used in this application:**