


Methadone by MSD on a ZB-1MS**Column:** Zebtron[™] ZB-1MS, GC Cap. Column 10 m x 0.18 mm x 0.18 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 μ m**Order No:** 7CD-G011-08**Oven Profile:** 150 °C to 320 °C @ 20 °C/min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 10:1 1 μ L @ 250°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 24ppm
Sim ions:
(Quant ion*) D9-Methadone (*78,303): Methadone (*72,294,223)**Products used in this application:**

16692

**ANALYTES:**

- 1 Methadone-d9
- 2 Methadone

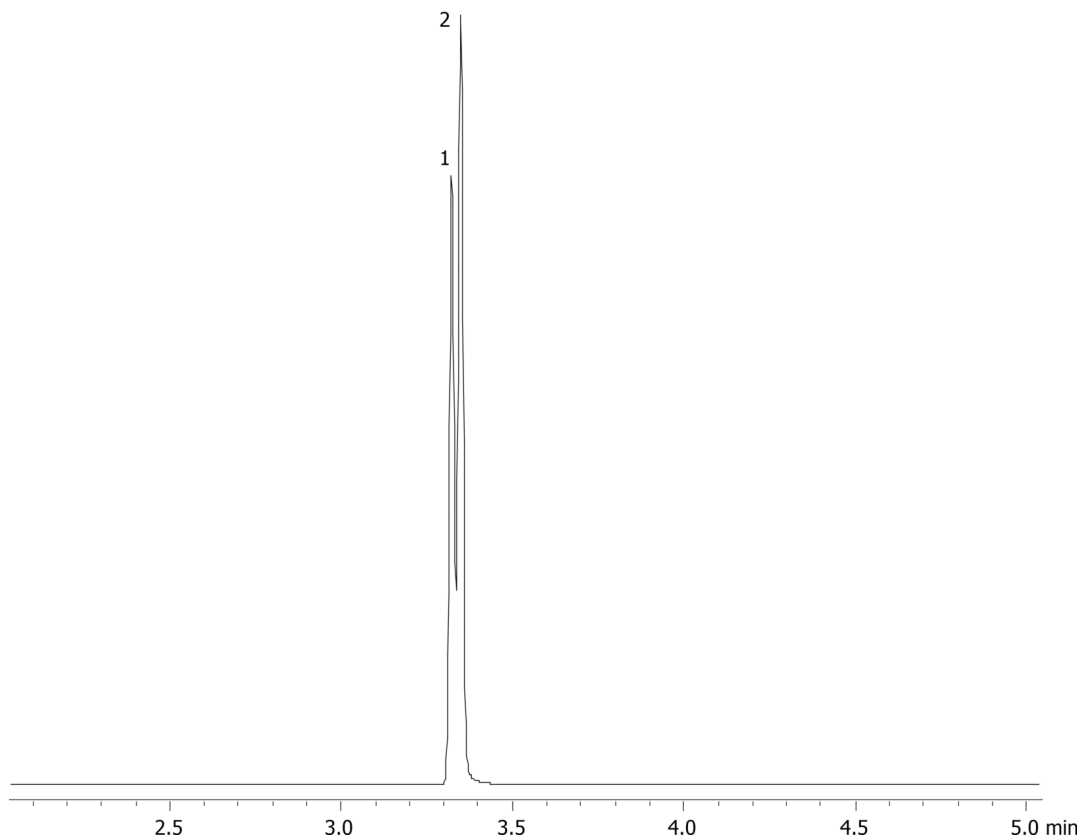