

Fluoxetine on Lux 3 μ m Cellulose-1**Column:** Lux[®] 3 μ m Cellulose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4458-E0**Elution Type:** Isocratic**Eluent A:** 0.1% DEA in Hexane**Eluent B:** 0.1% DEA in Isopropanol

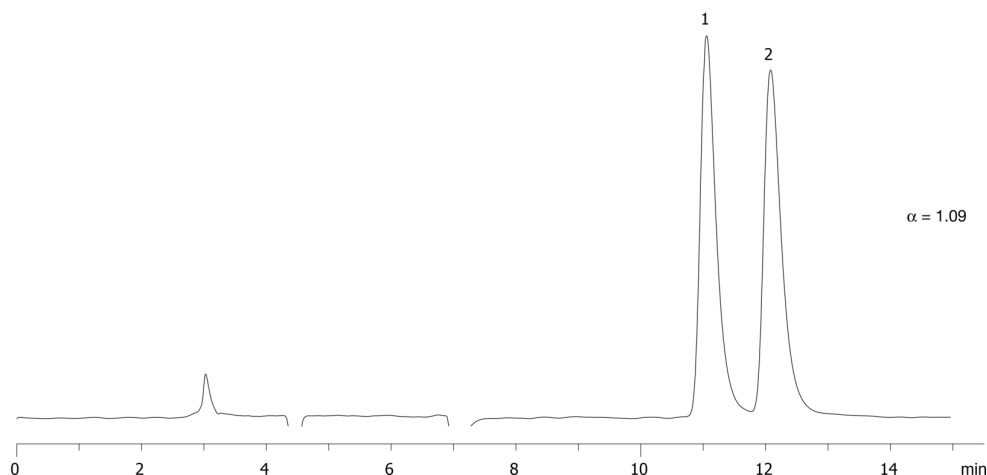
Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5

Flow Rate: 1 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 230 nm (ambient)**Detector Info:** α =1.12 and R_s =2.07**Analyst Note:** SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm, 10/Pk Part No.: AJ0-8403

- Holder Part No.: KJ0-4282

16731

**Products used in this application:****ANALYTES:**

- 1 Fluoxetine
- 2 Fluoxetine

