


**Sterols from Lard by FID on ZB-5HT****Column:** Zebtron<sup>™</sup> ZB-5HT, GC Cap. Column 30 m x 0.25 mm x 0.10 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.1 µm**Order No:** 7HG-G015-02**Oven Profile:** 220 °C to 350 °C @ 15 °C/min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Splitless :1 0.5 µL @ 350°C**Detection:** Flame Ionization (FID) (350°C)**Analyst Note:** Sample Prep: Saponification, Solid Phase Extraction (SPE), and Derivatization.  
Recommended Accessories:  
Recommended Liner: Zebtron PLUS Single Taper Z-Liner, 4 mm ID  
Liner Part No.: AG2-0A13-05 (for Agilent systems)**Products used in this application:**

16752

**ANALYTES:**

- 1 Cholesterol
- 2 Betulin

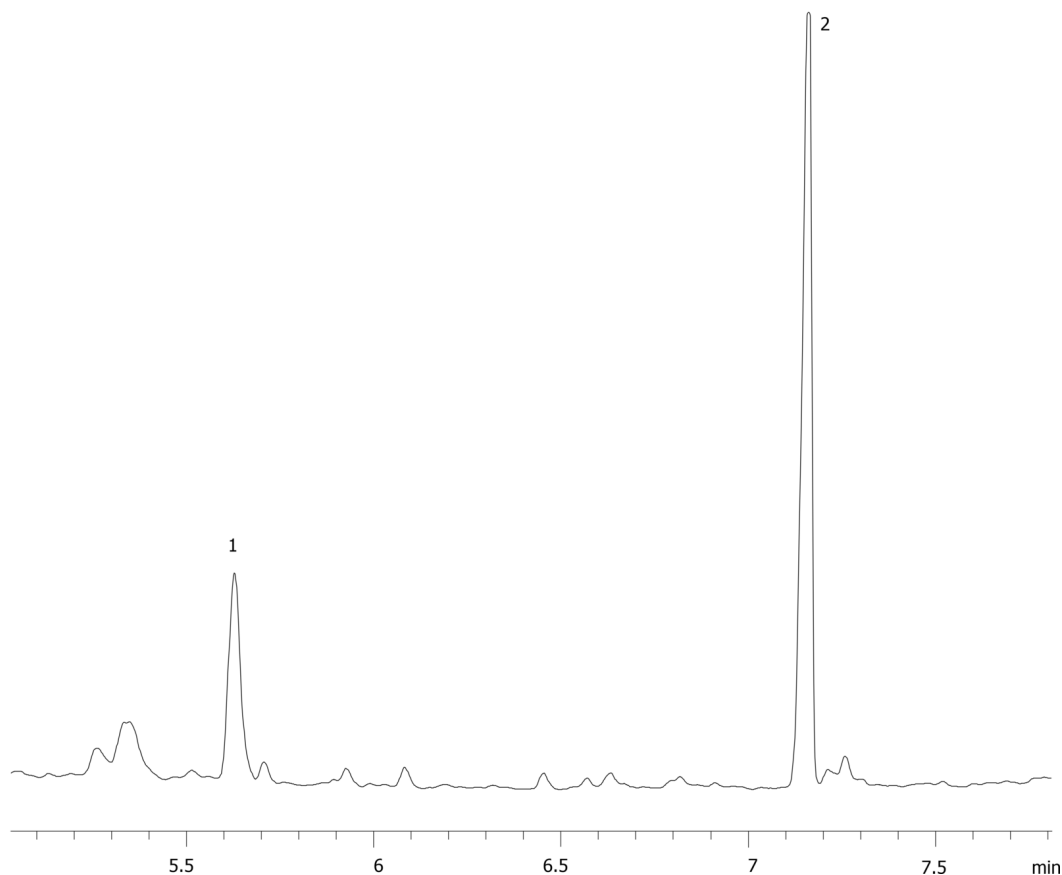