


Sterols from Lard by FID on ZB-MR2**Column:** Zebtron[™] ZB-MultiResidue[™]-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2 μ m**Order No:** 7HG-G017-10**Oven Profile:** 220 °C to 320 °C @ 13 °C/min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Splitless :1 0.5 μ L @ 300°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Sample Prep: Saponification, Solid Phase Extraction (SPE), and Derivatization.**Products used in this application:**

16757

**ANALYTES:**

- 1 Cholesterol
- 2 Stigmasterol
- 3 Betulin

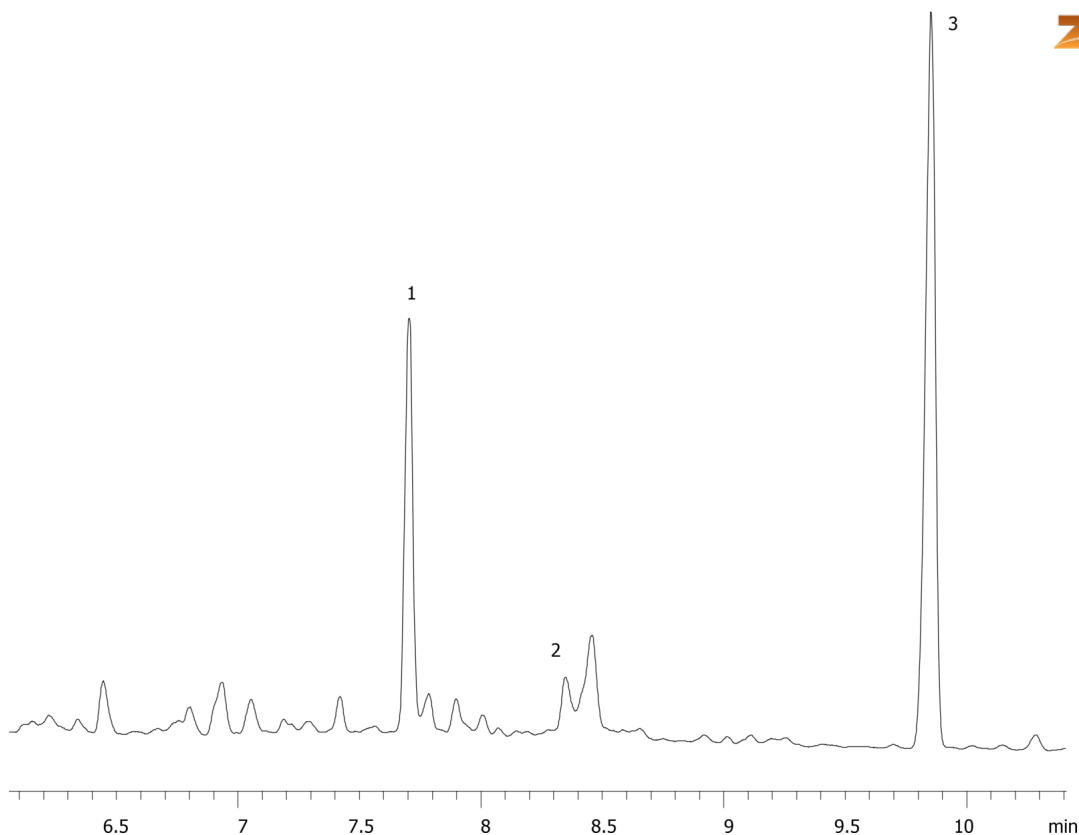