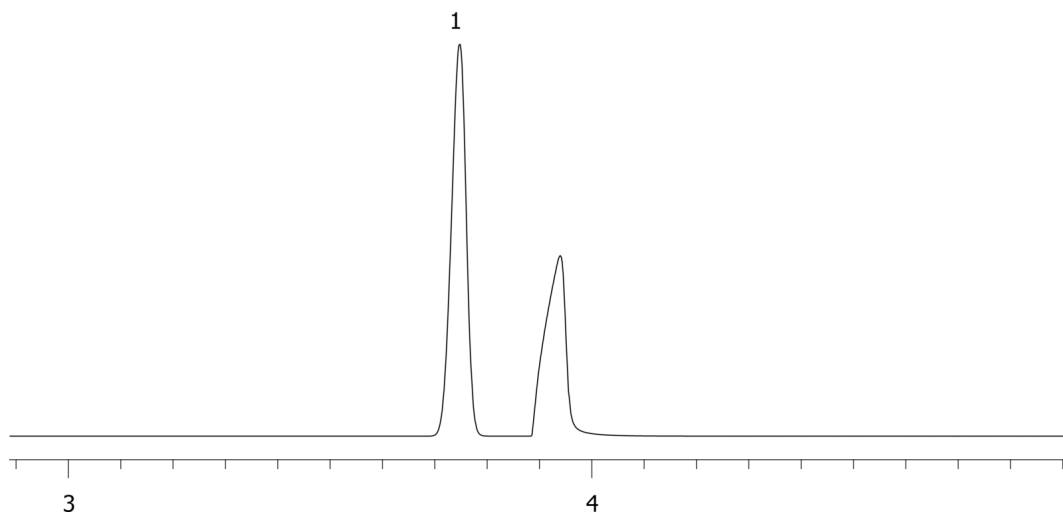


Benzene and butanol by MSD on a ZB-5MSi (a)**Column:** Zebtron[™] ZB-5MSi, GC Cap. Column 30 m x 0.25 mm x 0.50 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.5 μ m**Order No:** 7HG-G018-17**Oven Profile:** 30°C isothermal**Carrier Gas:** Constant Flow Helium, 55 cm/sec**Injection:** Split 50:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Analytes are 50ppm in DMSO
(Bad peak shape, not for view by sales)**Products used in this application:**

17684

**ANALYTES:****1** Benzene**2** Butanol