

Preparative Purification of Methocarbamol on Lux Cellulose 1

Column: Lux® 5 µm Cellulose-1, LC Column 100 x 21.2 mm, AXIA™ Packed, Ea

Dimensions: 100 x 21.2 mm ID

Order No: 00D-4459-P0-AX

Elution Type: Isocratic

Eluent A: Methanol

Eluent B: Isopropanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	90	10	

Flow Rate: 20 mL/min

Col. Temp.: ambient


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Detector Info: 32 mg

Analyst Note: Injection: 32 mg in 640 µL. cf: AppIDs 17780 (5 µg in 2 µL on 100 x 4.6 mm ID) and 17782 (80 mg in 1600 µL on 250 x 21.2 mm)



Products used in this application:

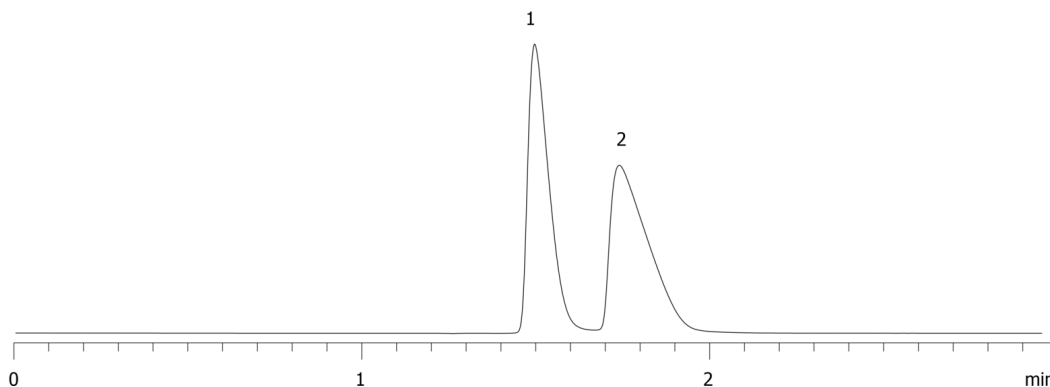

ANALYTES:

- 1 Methocarbamol
Retention Time: 1.5 min
- 2 Methocarbamol
Retention Time: 1.74 min

