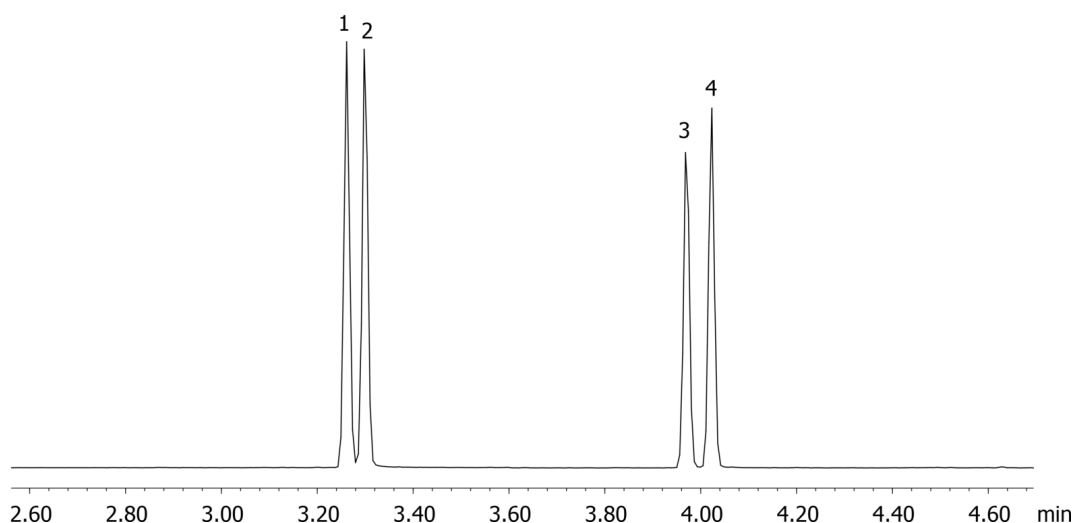


Amphetamines by MSD on a ZB-Drug1**Column:** Zebtron[™] ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:****Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:** 110°C to 320°C @ 25 °C/min**Carrier Gas:** Constant Flow Helium, 55 cm/sec**Injection:** Split 10:1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 10ppm derivatized with HFBA

Sim ions:

(Quant ion*) D11-Amphetamine (*244,128): Amphetamine (*240,118,117): D14-Methamphetamine (*261,213):
Methamphetamine (*254,210,118)**Products used in this application:**

17807

**ANALYTES:**

- 1 Amphetamine-d11
- 2 Amphetamine
- 3 Methamphetamine-d14
- 4 Methamphetamine

