


**Blood Alcohols on ZB-BAC-1****Column:** Zebtron<sup>™</sup> ZB-BAC-1, GC Cap. Column 30 m x 0.32 mm x 1.80  $\mu$ m, Ea**Phase:** Proprietary Blood Alcohol Phase**Dimensions:** 30 meters x 0.32 mm x 1.8  $\mu$ m**Order No:** 7HM-G021-31**Oven Profile:** 40 °C Isothermal**Carrier Gas:** Constant Flow Helium, 80 cm/sec**Injection:** Split 5:1 0.5  $\mu$ L @ 140°C**Detection:** Flame Ionization (FID) (220°C)**Analyst Note:** 0.050% Blood Alcohols Components in Water (1mL in 20mL Headspace Vial)

17897

**Products used in this application:****ANALYTES:**

- 1 Methanol
- 2 Acetaldehyde
- 3 Ethanol
- 4 Isopropanol
- 5 Acetone
- 6 t-Butanol
- 7 2-Butanol

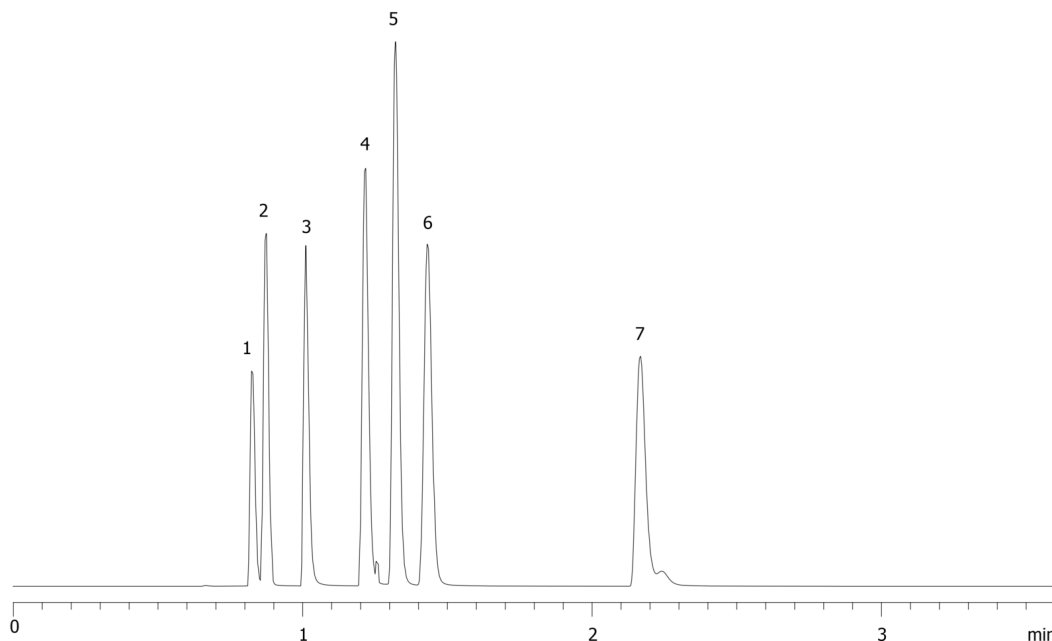