

EPA Method 8330: Explosive Mix Using Luna 5u PFP

Column: Luna® 5 µm PFP(2) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4448-E0

Elution Type: Gradient

Eluent A: ACN

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	30	70
	2	2	30	70
	3	22	50	50
	4	22.01	30	70
	5	25	30	70

Flow Rate: 1


Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, PFP 4 x 3.0mm, 10/Pk Part No.: AJ0-8327

- Holder Part No.: KJ0-4282

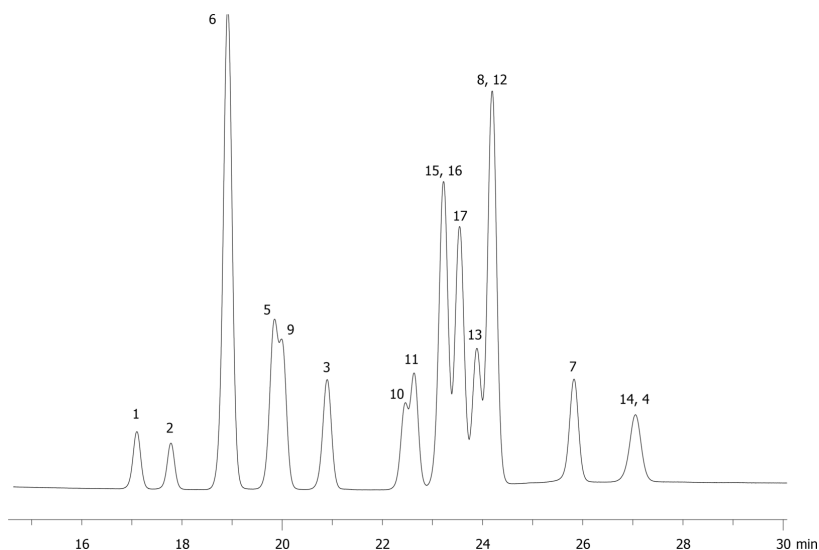

Products used in this application:



EPA Method 8330: Explosive Mix Using Luna 5u PFP

ANALYTES:

- 1** Octahydro-1,3,5,7-tetranitro-1,3,5,7-tetrazocine (HMX)
- 2** Hexahydro-1,3,5-trinitro-1,3,5-triazine
- 3** 1,3,5-Trinitrobenzene
- 4** Methyl-2,4,6-trinitrophenylnitramine (Tetryl)
- 5** 1,3-Dinitrobenzene
- 6** Nitrobenzene
- 7** Nitroglycerin
- 8** 2,4,6-Trinitrotoluene
- 9** 3,5-Dinitroaniline
- 10** 4-Amino-2,6 dinitrotoluene
- 11** 2 Amino-4,6 dinitrotoluene
- 12** 2,6-Dinitrotoluene
- 13** 2,4-Dinitrotoluene
- 14** Pentaerythritol tetranitrate
- 15** 2-Nitrotoluene
- 16** 4-Nitrotoluene
- 17** 3-Nitrotoluene

