

Narcotics on Kinetex 2.6µm C18 50 x 4.6mm ID

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4462-E0

Elution Type: Gradient

Eluent A: 10mM NH4OAc pH5.5

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.34	95	5
	3	2.61	60	40


Flow Rate: 2.94 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
18100



Products used in this application:

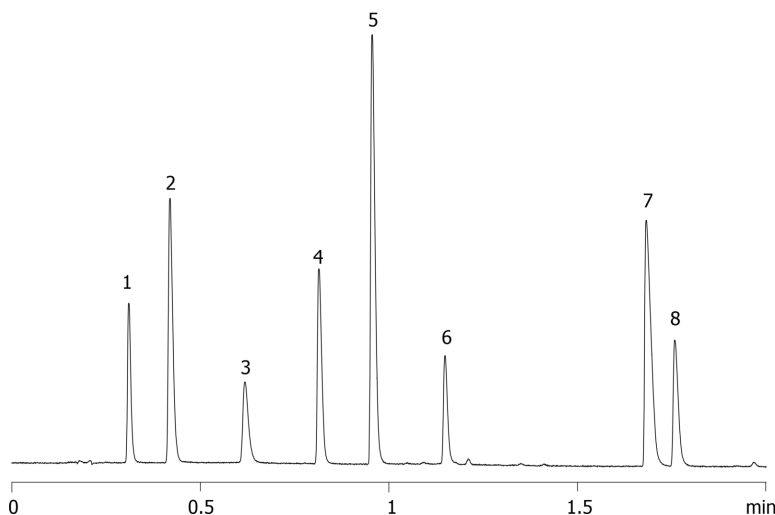

ANALYTES:

- 1 Normorphine
Retention Time: 0.311 min
- 2 Morphine
Retention Time: 0.421 min
- 3 Hydromorphone
Retention Time: 0.619 min
- 4 Norcodeine
Retention Time: 0.816 min
- 5 Codeine
Retention Time: 0.967 min
- 6 Hydrocodone
Retention Time: 1.15 min
- 7 Cocaine
Retention Time: 1.685 min
- 8 Norcocaine
Retention Time: 1.76 min

