

## Pharma ABN Test 12 Mix on Luna 5um C18(2) 250x4.6mm 1mL/min

**Column:** Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4252-E0

**Elution Type:** Gradient

**Eluent A:** water w/ 0.1% formic acid

**Eluent B:** acetonitrile w/ 0.1% formic acid

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 95    | 5     |
|                   | 2        | 2.79       | 95    | 5     |
|                   | 3        | 36.15      | 5     | 95    |
|                   | 4        | 36.38      | 95    | 5     |
|                   | 5        | 50.93      | 95    | 5     |

**Flow Rate:** 1 mL/min

**Col. Temp.:** 45 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (25 °C)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287


- Holder Part No.: KJ0-4282



Products used in this application:



18230

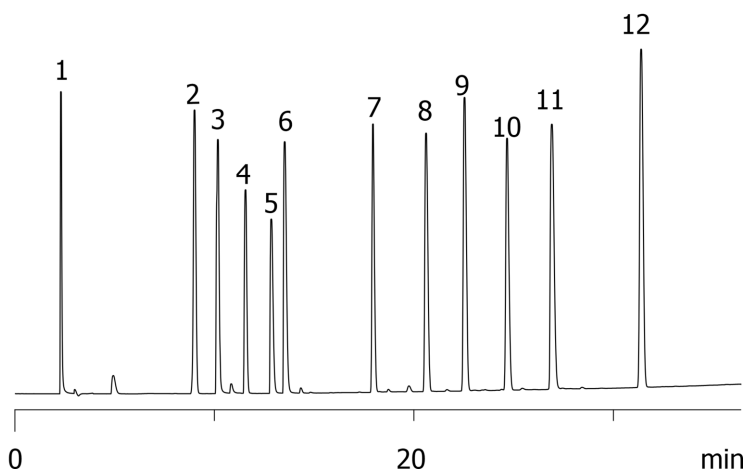

## Pharma ABN Test 12 Mix on Luna 5um C18(2) 250x4.6mm 1mL/min

### ANALYTES:

- 1** Pyridine
- 2** Acetaminophen
- 3** Quinine
- 4** Acebutolol
- 5** Chlorpheniramine
- 6** Triprolidine
- 7** Prednisolone
- 8** 4-Chlorobenzoic acid
- 9** 4-Chlorocinnamic acid
- 10** Diazepam
- 11** Diflunisal
- 12** Hexanophenone

