

## Agilent 1200 - Kinetex 2.6µm C18 50 x 4.6mm ID (2.5 mL/min)

**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 4.6 mm, Ea

**Dimensions:** 50 x 4.6 mm ID

**Order No:** 00B-4462-E0

**Elution Type:** Isocratic

**Eluent A:** water

**Eluent B:** acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 50    | 50    |


**Flow Rate:** 2.5 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



Products used in this application:

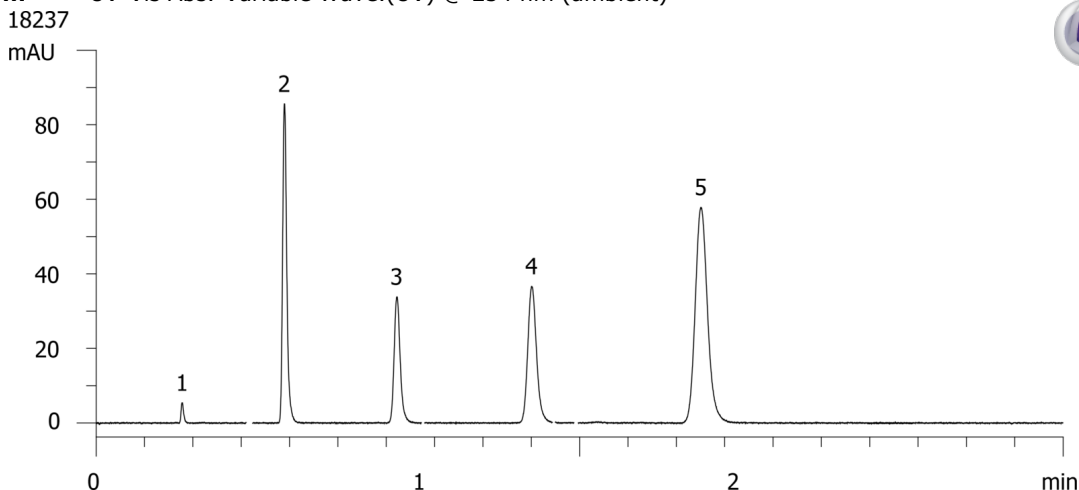

### ANALYTES:

- 1** Uracil  
Retention Time: 0.178 min
- 2** Acetophenone  
Retention Time: 0.39 min
- 3** Benzene  
Retention Time: 0.622 min
- 4** Toluene  
Retention Time: 0.901 min
- 5** Naphthalene  
Retention Time: 1.251 min

