

Extraction of Sulfa Drugs from Porcine Plasma using Strata X-C

Column: Synergi™ 4 µm Max-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4337-E0

Elution Type: Gradient

Eluent A: 20 mM KH₂PO₄ (pH2.5)

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	15	20	80

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282



18293

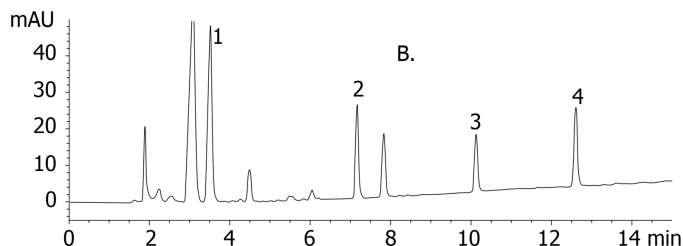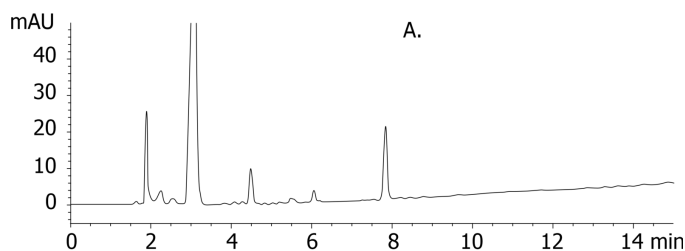

ANALYTES:

1 Sulfanilamide

Retention Time: 3.49 min

2 Sulfathiazole

Retention Time: 7.14 min

3 Sulfamethoxazole

Retention Time: 10.11 min

4 Sulfaquinoxaline

Retention Time: 12.58 min

