

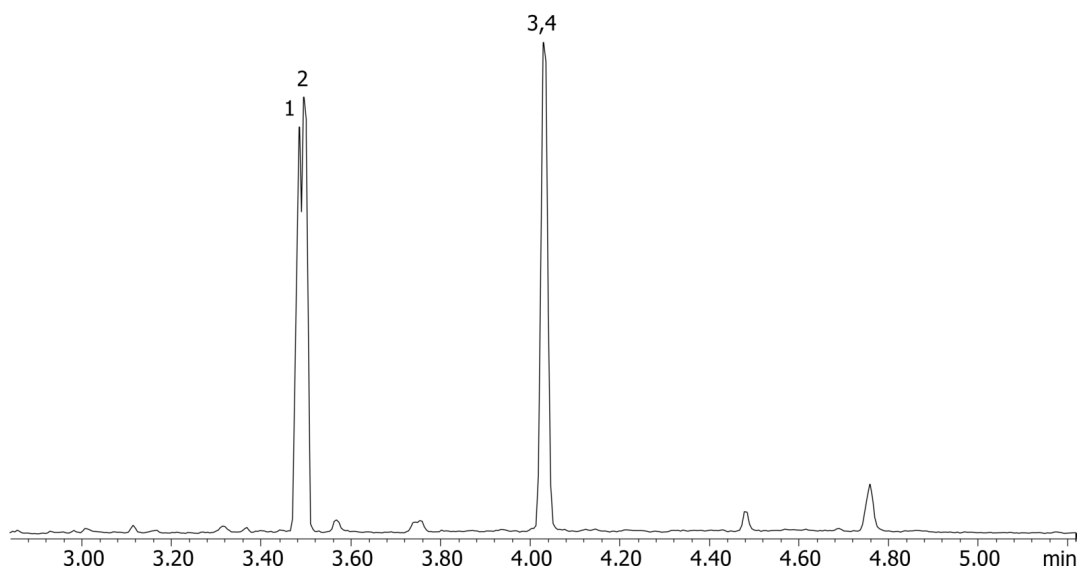
Opiates from urine by Strata-XC on a ZB-Drug1**Column:** Zebtron™ ZB-Drug-1, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea**Phase:** Proprietary Drug Phase**Dimensions:** 10 meters x 0.18 mm x 0.18 µm**Order No:** 7CD-G023-08**Oven Profile:** 185 °C to 320 °C @ 30 °C/min for 3 min**Carrier Gas:** Constant Flow Helium, 0.7 mL/min**Injection:** Split 10:1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Analytes are 300 ng/mL in urine derivatized with propionic anhydride

Sim ions:

(Quant ion*) D6-Codeine (*361,288,235): Codeine (*355,282,229): D3-Morphine (*344,271,400): Morphine (*341,268,397)

**Products used in this application:**

18301

**ANALYTES:**

- 1 Codeine-d6
- 2 Codeine
- 3 Morphine-d3
- 4 Morphine

