

SecurityGuard Comparisons on BioSep SEC 3000 300 x 7.8mm Without Cartridge

Column: BioSep™ 5 µm SEC-s3000 290 Å, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-2146-K0

Elution Type: Isocratic

Eluent A: 50mM Phosphate + 300mM NaCl pH=6.8 Buffer


Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min


Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: Analysis performed without SecurityGuard Cartridge (cf. App IDs: 18793 and 18794) using Aqueous SEC Test Mix (PN: A10-3042)
18792

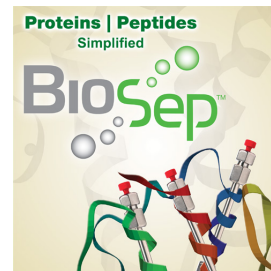

Products used in this application:

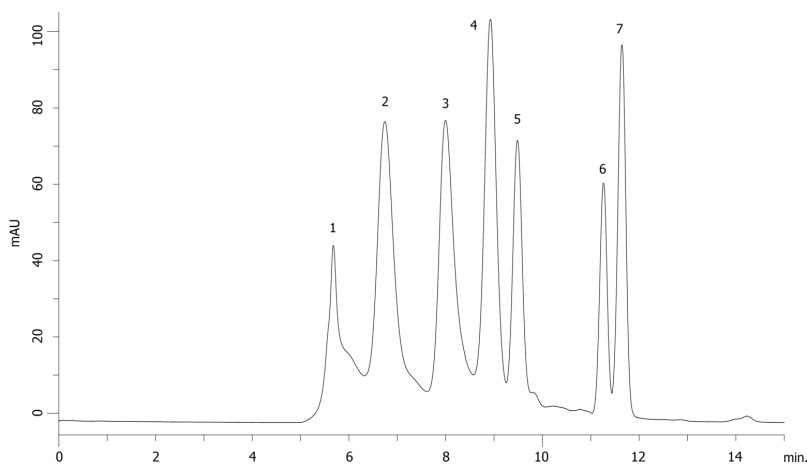

ANALYTES:

- 1 IgM Impurity
- 2 IgM
Retention Time: 6.74 min
- 3 Thyroglobulin
Retention Time: 7.89 min
- 4 Beta Amylase
Retention Time: 8.92 min
- 5 Ovalbumin
Retention Time: 9.48 min
- 6 Myoglobin
Retention Time: 11.26 min
- 7 Uridine
Retention Time: 11.64 min

