

## Tryptic digest of BSA using Kinetex 2.6u C18 150 x 4.6mm ID

**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4462-E0

**Elution Type:** Gradient

**Eluent A:** Water with 0.1% TFA

**Eluent B:** Acetonitrile with 0.1% TFA

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 97    | 3     |
|                   | 2        | 30         | 50    | 50    |

**Flow Rate:** 1.2 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)

19078



Products used in this application:

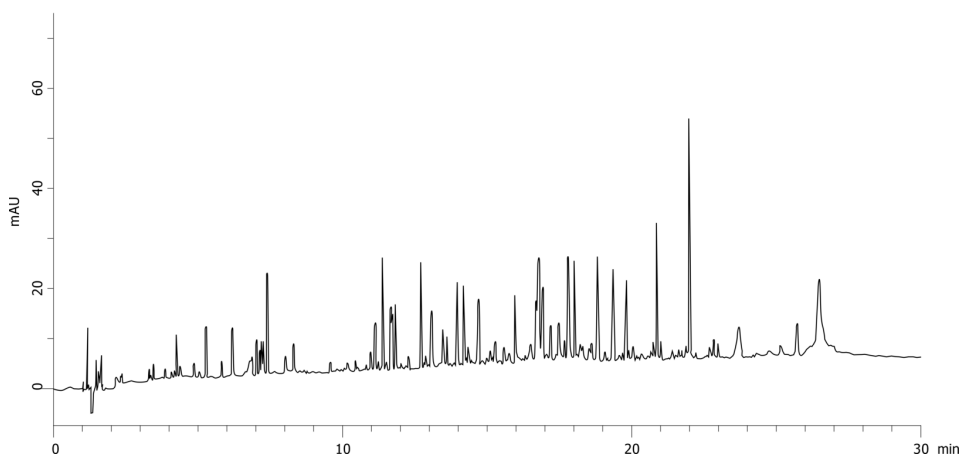

### ANALYTES:

- 1 Bovine Serum Albumin (BSA)

