

## CD Detection of TSO on LUX 5µm Cellulose-1 250 x 4.6 mm ID in NP

**Column:** Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4459-E0

**Elution Type:** Isocratic

**Eluent A:** 0.1% DEA in Hexane

**Eluent B:** 0.1% DEA in IPA

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 10         | 90    | 10    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** 25 °C

**Detection:** Circular Dichroism @ 254 nm (25 °C)


**Detector Info:** Jasco CDD

**Analyst Note:** Please see the corresponding CD trace (attached document).

19149



Products used in this application:

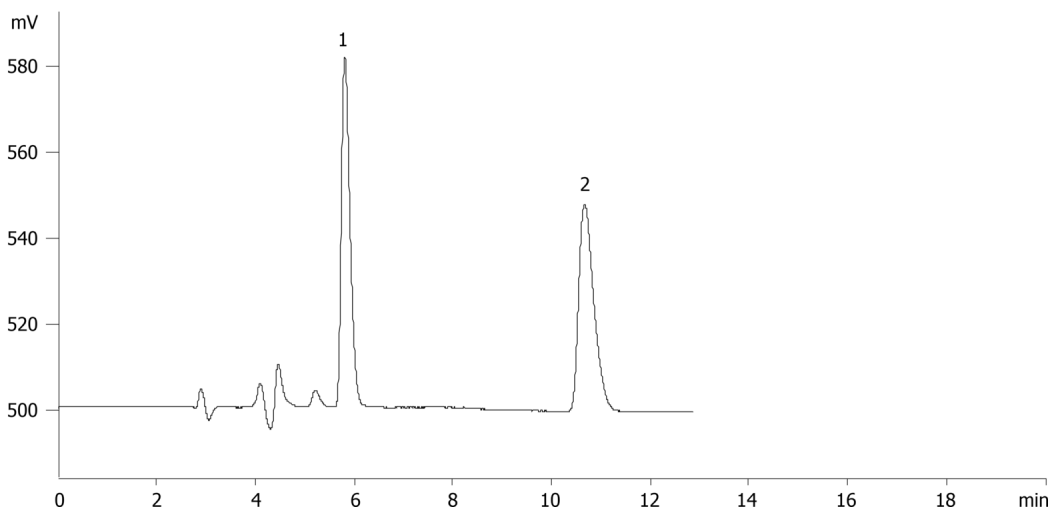

### ANALYTES:

- 1 trans-Stilbene oxide
- 2 trans-Stilbene oxide

