HPLC Application

ID No.: **19167**



Pain Panel with and without sample focusing on Kinetex C18 2.6µm 50 x 2.1 mm by LC/MS

Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea


50 x 2.1 mm ID **Dimensions:** Order No: 00B-4462-AN Elution Type: Gradient

Eluent A: 0.1% Formic acid in H2O Eluent B: 0.1% Formic acid ACN

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	0.5	95	5
	3	2.5	5	95
	4	3.5	5	95
	5	3.51	95	5
	6	5.3	95	5

Flow Rate: 450 µL/min Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ amu (ambient)

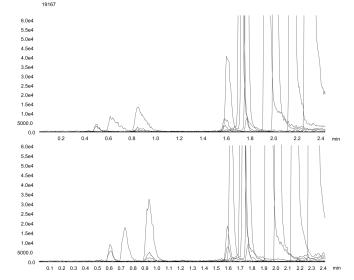



Products used in this application:





Page 1 of 2

HPLC ApplicationID No.: 19167



Pain Panel with and without sample focusing on Kinetex C18 2.6 μ m 50 x 2.1 mm by LC/MS

ANALYTES:

- Morphine 1
- 2 Oxymorphone
- 3 Hydromorphone
- Codeine 4
- 5 Naloxone
- 6 Oxycodone
- 7 Naltrexone
- 8 Hydrocodone
- 9 Tramadol
- 10 Meperidine
- 11 Fentanyl
- 12 Buprenorphine
- **13** Methadone