


Narcotics 8-mix at Neutral pH on Kinetex 2.6 µm C18, 50 x 2.1mm**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4462-AN**Elution Type:** Gradient**Eluent A:** 10mM NH₄OAc pH5.5**Eluent B:** ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.25	95	5
	3	2	60	40
	4	2.01	95	5
	5	3	95	5

Flow Rate: 0.8 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Detector Info:** 550 bar back pressure**Products used in this application:****ANALYTES:**

- 1 Normorphine
- 2 Morphine
- 3 Hydromorphone
- 4 Norcodeine
- 5 Codeine
- 6 Hydrocodone
- 7 Cocaine
- 8 Norcocaine

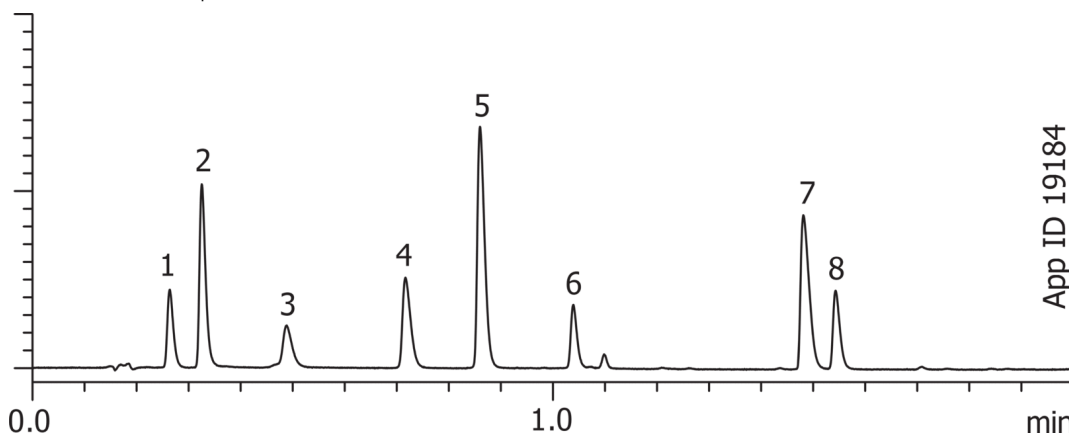