

## Insulin impurity profiling using Kinetex 2.6µm C18 150x4.6mm

**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4462-E0

**Elution Type:** Gradient

**Eluent A:** 20mM Potassium phosphate pH 2.7

**Eluent B:** Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 78    | 22    |
|                   | 2        | 30         | 66    | 34    |

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)

19277



Products used in this application:

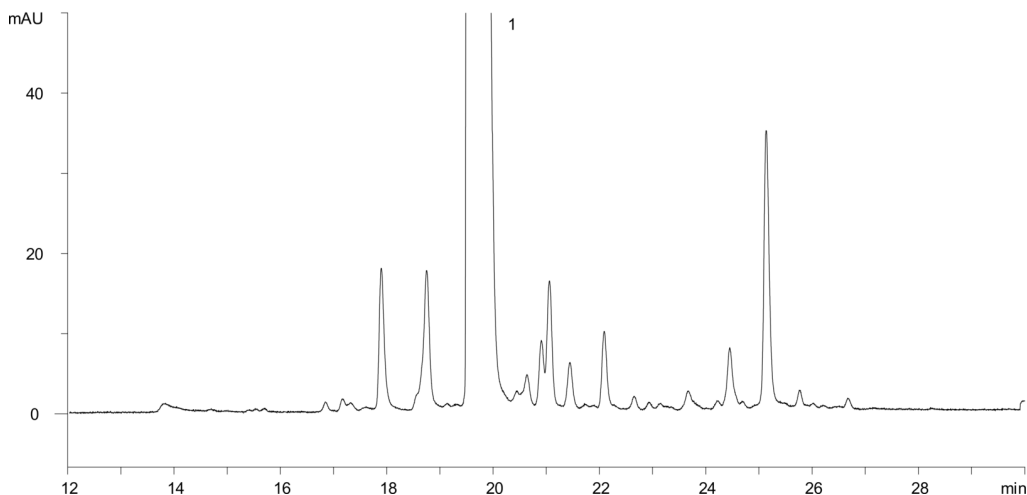

### ANALYTES:

1 Insulin

