

EPA Method 8151 - Chlorinated Herbicides on a ZB-35-HT 30x0.32x0.25 (Dual column)

Column: Zebtron™ ZB-35HT, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea
Phase: 35% Phenyl 65% Dimethylpolysiloxane
Dimensions: 30 meters x 0.32 mm x 0.25 µm
Order No: 7HM-G025-11
Oven Profile:

Carrier Gas: Constant Flow Helium, 2.4 mL/min
Injection: Pulsed 0.000000000:1 0 µL @ 250°C
Detection: Refractive Index (350°C)
Analyst Note: Analytes are 100ppb in isooctane.


Dual column with ZB-XLB-HT



Products used in this application:



19282

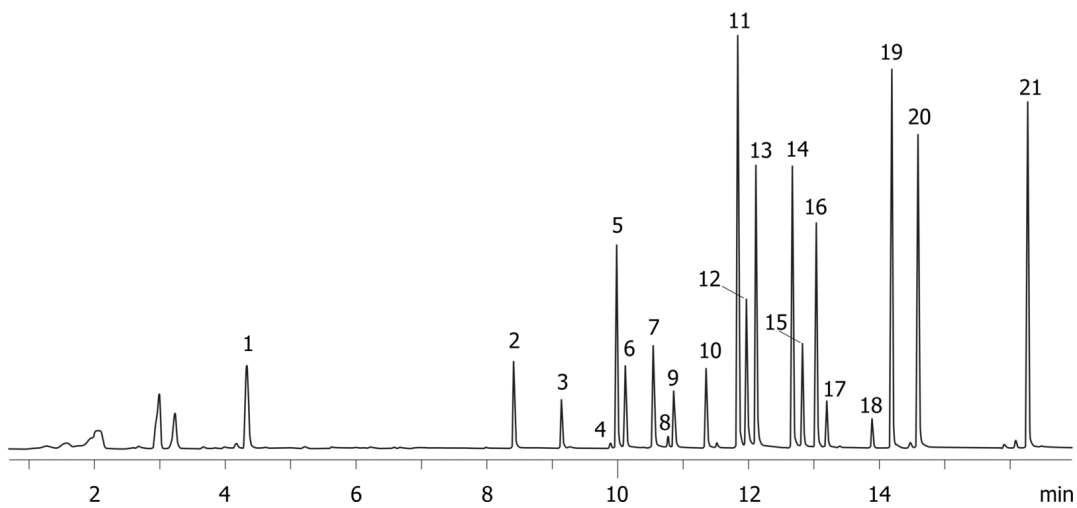

EPA Method 8151 - Chlorinated Herbicides on a ZB-35-HT 30x0.32x0.25 (Dual column)

ANALYTES:

- 1 Dalapon
- 2 3,5-Dichlorobenzoic acid
- 3 4-Nitrophenol
- 4 DCAA
- 5 Dicamba
- 6 MCPP
- 7 MCPA
- 8 DBOB
- 9 Dichlorprop
- 10 2,4-D
- 11 Pentachlorophenol
- 12 contaminant
- 13 Silvex
- 14 Chloramben
- 15 2,4,5-T
- 16 Dinoseb
- 17 2,4-DB
- 18 Bentazon
- 19 Picloram
- 20 Dacthal
- 21 Acifluorfen

