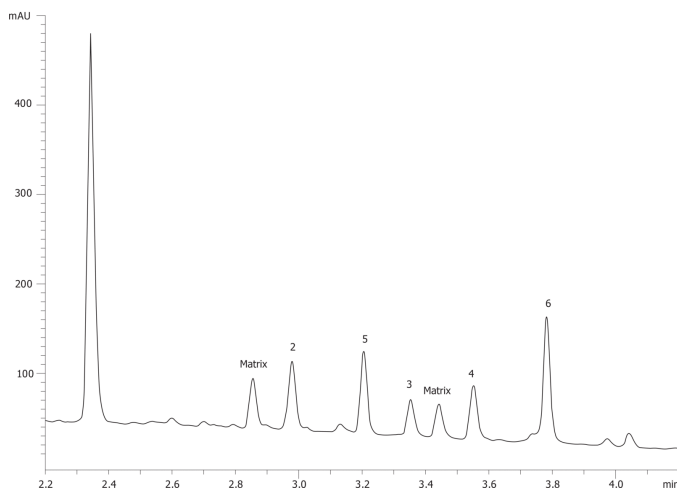


Veterinary Drugs in Horse Urine on Kinetex 2.6µm C8 50x4.6mm @ 1.50mL/min**Column:** Kinetex® 2.6 µm C8 100 Å, LC Column 50 x 4.6 mm, Ea**Dimensions:** 50 x 4.6 mm ID**Order No:** 00B-4497-E0**Elution Type:** Gradient**Eluent A:** 0.1% Phosphoric Acid in Water**Eluent B:** 0.1% Phosphoric Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	2	50	50
	3	3.75	30	70
	4	4.5	5	95
	5	5	70	30

Flow Rate: 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)

19472

**Products used in this application:****ANALYTES:**

- 1 Naproxen
- 2 Fenoprofen
- 3 Ibuprofen
- 4 Phenylbutazone
- 5 Indomethacin
- 6 Meclofenamic acid
- 7 matrix impurity
- 8 matrix impurity

