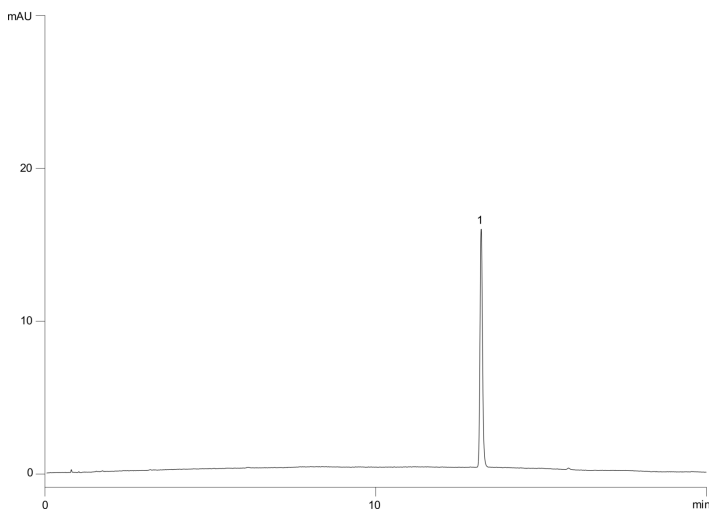


USP Hydrocortisone Acetate on Kinetex 2.6µm XB-C18 75x4.6mm**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 75 x 4.6 mm, Ea**Dimensions:** 75 x 4.6 mm ID**Order No:** 00C-4496-E0**Elution Type:** Gradient**Eluent A:** 80:20 Water:Acetonitrile**Eluent B:** 70:30 Acetonitrile:Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	5	90	10
	3	25	10	90
	4	30	10	90
	5	35	90	10
	6	40	90	10

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19517

**ANALYTES:**

- 1 Hydrocortisone acetate

**Products used in this application:**