

# HPLC Application

ID No.: 19562

## USP Loratadine on Kinetex 2.6µm C8 150x4.60mm

**Column:** Kinetex® 2.6 µm C8 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4497-E0

**Elution Type:** Isocratic

**Eluent A:** 31.5:31.5:36.8 Methanol:Acetonitrile:Solution A(1.74mg/mL Anhydrous Dibasic Potassium)

**Gradient**

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 15         | 100   |


**Profile:**

**Flow Rate:** 1

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

19562

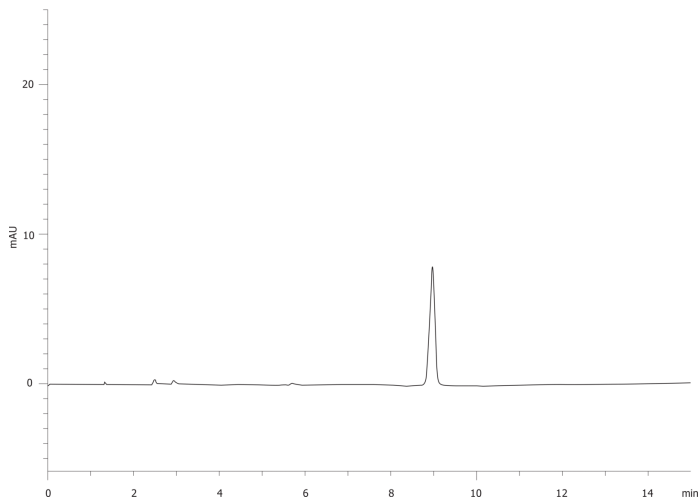

### ANALYTES:

1 Loratadine



Products used in this application:

