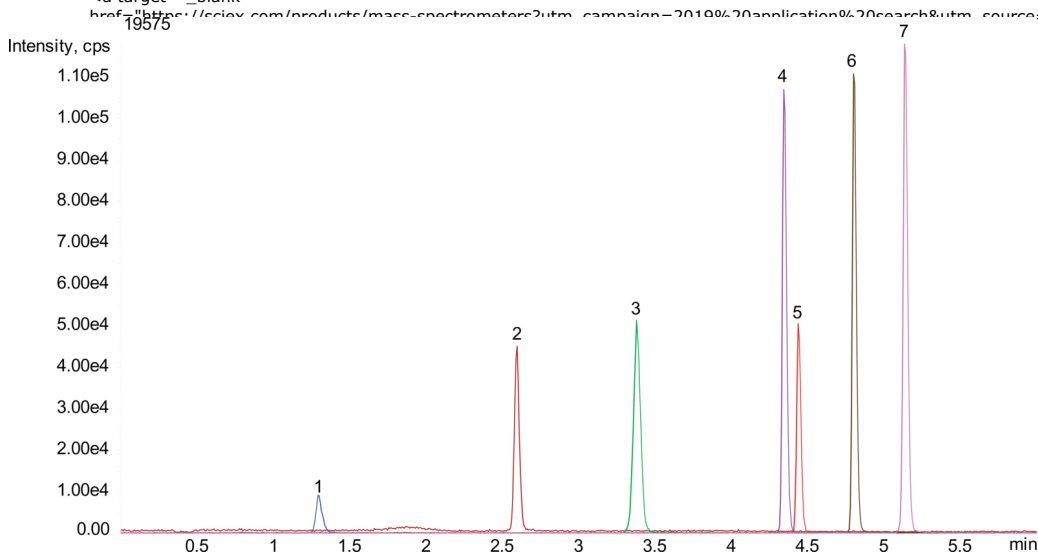


Triazine Pesticides (EPA 536) on Kinetex 2.6µm XB-C18 by LC/MS/MS**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4496-AN**Elution Type:** Gradient**Eluent A:** 5 mM Ammonium Acetate**Eluent B:** 100% MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.25	60	40
	3	2	60	40
	4	3	25	75
	5	4	25	75
	6	4.1	95	5

Flow Rate: 0.3 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Detector Info:** SCIEX<**ANALYTES:****1** Atrazine-desethyl-desisopropyl

Retention Time: 1.29 min

2 Atrazine-desisopropyl

Retention Time: 2.59 min

3 Atrazine-desethyl

Retention Time: 3.38 min

4 Cyanazine

Retention Time: 4.38 min

5 Simazine

Retention Time: 4.44 min

6 Atrazine

Retention Time: 4.8 min

7 Propazine

Retention Time: 5.14 min

**Products used in this application:**