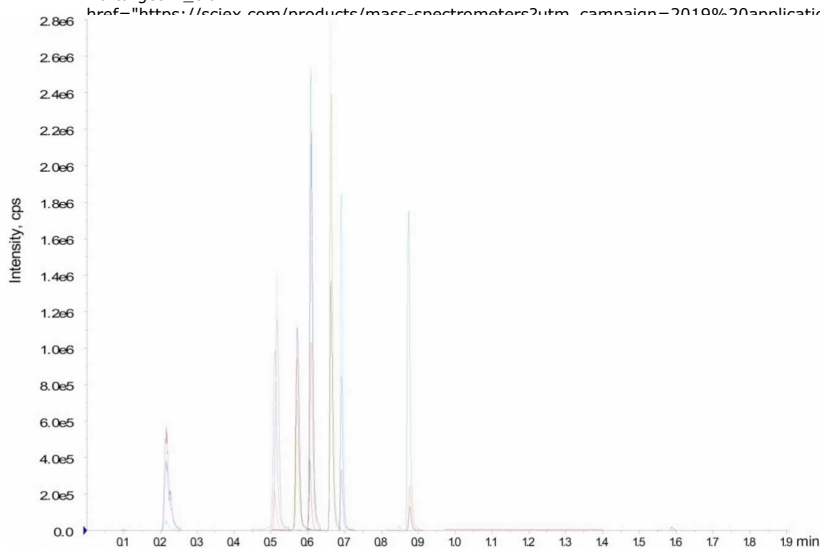


Kinetex SAMHSA LC/MS Run**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4462-AN**Elution Type:** Gradient**Eluent A:** 0.1% Formic Acid**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1.7	5	95
	3	1.9	5	95
	4	1.91	95	5

Flow Rate: 1.2 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Detector Info:** <a target="_blank"**Products used in this application:**

Kinetex SAMHSA LC/MS Run

ANALYTES:

- 1** Morphine
- 2** Morphine-d3
- 3** Codeine
- 4** Codeine-d6
- 5** Amphetamine
- 6** Methamphetamine
- 7** MDA
- 8** MDA-d5
- 9** 6-MAM
- 10** MAM-d3
- 11** MDMA
- 12** MDMA-d5
- 13** MDEA
- 14** MDEA-d5
- 15** Benzoylecgonine
- 16** Benzoylecgonine-d3
- 17** PCP
- 18** PCP-d5

