

## Clopidogrel Bisulfate on Lux 5µm Cellulose-3 in PO

**Column:** Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4493-E0

**Elution Type:** Isocratic

**Eluent A:** MeOH (100%) DEA (0.1%)

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Lux Cellulose-3 4 x 3.0 mm, 10/Pk Part No.: AJ0-8622

- Holder Part No.: KJ0-4282

19646



Products used in this application:

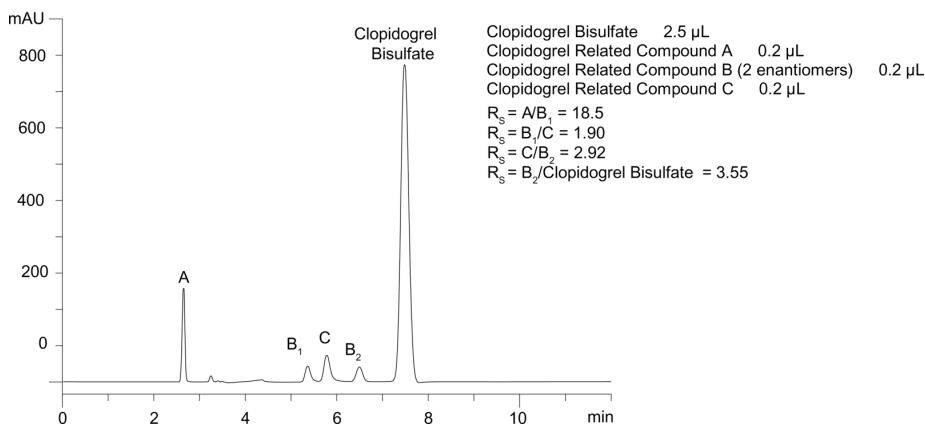

### ANALYTES:

- 1 Clopidogrel related compound A  
Retention Time: 2.66 min
- 2 Clopidogrel related compound B1  
Retention Time: 5.363 min
- 3 Clopidogrel related compound C  
Retention Time: 5.781 min
- 4 Clopidogrel related compound B2  
Retention Time: 6.494 min
- 5 Clopidogrel Bisulfate  
Retention Time: 7.5 min

