

## Steviol Glycoside standards on Kinetex 2.6u XB-C18 150x4.6mm

**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4496-E0

**Elution Type:** Isocratic

**Eluent A:** Water with 0.1% TFA

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	68	32

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** 70 °C

**Detection:** Circular Dichroism @ 0.0000000000 (ambient)

**Analyst Note:** SecurityGuard™ ULTRA Guard Cartridge System extends column lifetime.

- SecurityGuard ULTRA Cartridges, UHPLC C18 for 4.6mm ID Columns, 3/Pk Part No.: AJ0-8768


- Holder Part No.: AJ0-9000



Products used in this application:



19797

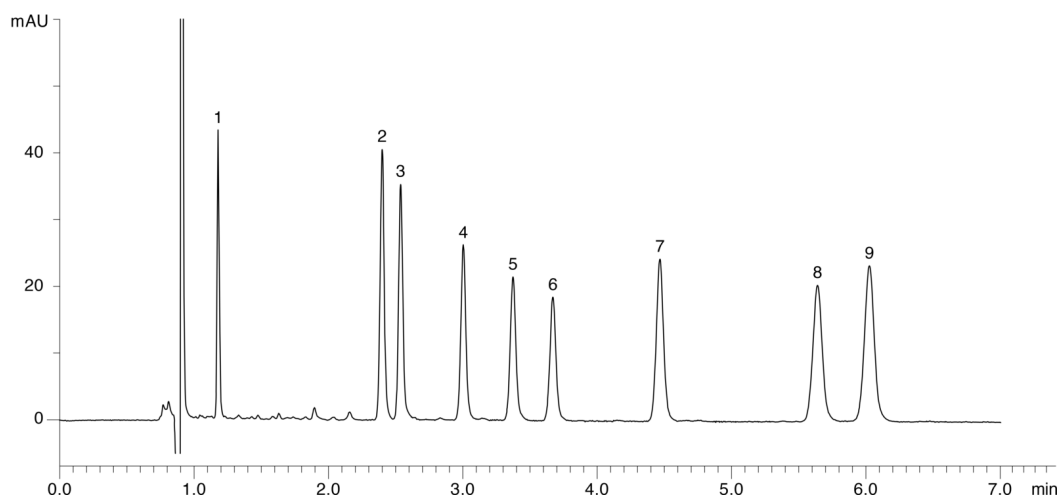

### ANALYTES:

- 1 Rebaudioside D
- 2 Rebaudioside A
- 3 Stevioside
- 4 Rebaudioside F
- 5 Rebaudioside C
- 6 Dulcoside A
- 7 Rubusoside
- 8 Rebaudioside B
- 9 Steviolbioside

