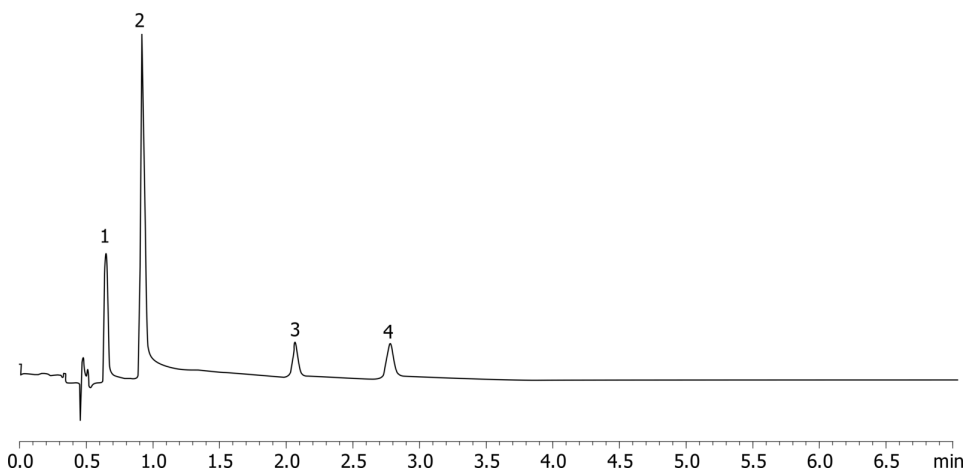


**USP Acetaminophen Acetylsalicylic Acid Caffeine on Kinetex 1.7u C18 100x2.10mm on Shimadzu Nexera****Column:** Kinetex® 1.7 µm C18 100 Å, LC Column 100 x 2.1 mm, Ea**Dimensions:** 100 x 2.1 mm ID**Order No:** 00D-4475-AN**Elution Type:** Isocratic**Eluent A:** 69:28:3 Water:Methanol:Acetic Acid**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 12 100**Flow Rate:** 0.4 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 275 nm (ambient)

19959

**Products used in this application:****ANALYTES:****1** Acetaminophen

Retention Time: 0.647 min

**2** Caffeine

Retention Time: 0.922 min

**3** Acetylsalicylic acid

Retention Time: 2.068 min

**4** Benzoic acid

Retention Time: 2.781 min

