


Yarra 3µm SEC-2000 EPO Aggregate Analysis

Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4512-K0
Elution Type: Isocratic
Eluent A: 50 mM Sodium Phosphate pH 6.8 / 0.3 M Sodium Chloride
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Circular Dichroism @ 0.0000000000 (ambient)

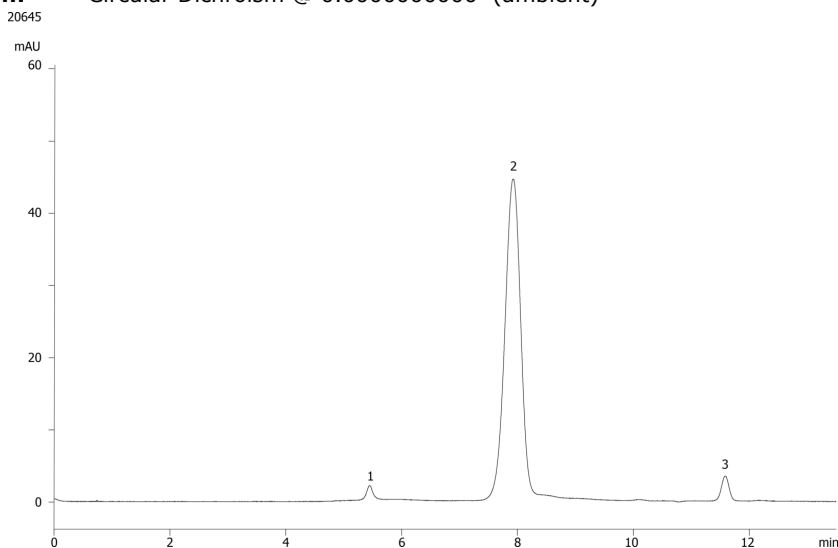

ANALYTES:

- 1** High MW aggregate
Retention Time: 5.4 min
- 2** EPO Monomer
Retention Time: 8 min
- 3** Low MW peak
Retention Time: 11.7 min



Products used in this application:

