

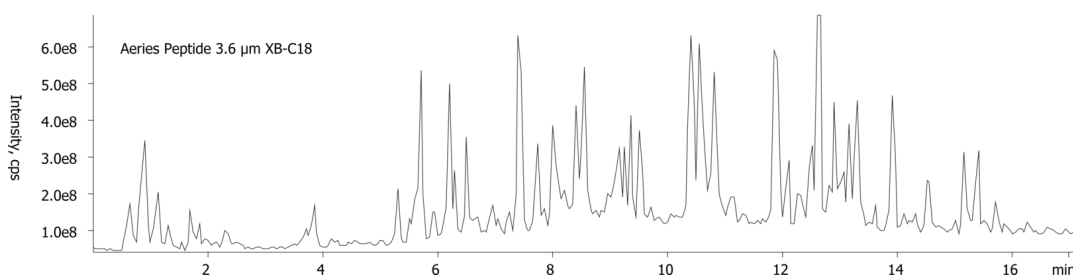
LC/MS Analysis of k-Casein Tryptic Digest with Aeris PEPTIDE 3.6u XB-C18**Column:** AerisTM 3.6 μ m PEPTIDE XB-C18 100 Å, LC Column 150 x 2.1 mm, Ea**Dimensions:** 150 x 2.1 mm ID**Order No:** 00F-4507-AN**Elution Type:** Gradient**Eluent A:** Water 0.1% Formic Acid**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	19	50	50
	3	21	5	95

Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Detector Info:** <a target="_blank"

20733

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019Q4%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

**Products used in this application:****ANALYTES:****1** k-Casein

Retention Time: 1 min

