


**Kinetex 5um - Mupirocin degradants using Kinetex 5um C18 250x4.6mm****Column:** Kinetex® 5µm C18 100 Å, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4601-E0**Elution Type:** Gradient**Eluent A:** Water with 0.1% TFA**Eluent B:** Acetonitrile with 0.1% TFA

| Gradient        | Step No. | Time (min) | Pct A | Pct B |
|-----------------|----------|------------|-------|-------|
| <b>Profile:</b> | <b>1</b> | 0          | 90    | 10    |
|                 | <b>2</b> | 40         | 30    | 70    |

**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**Products used in this application:****ANALYTES:**

- 1 Mupirocin degradants

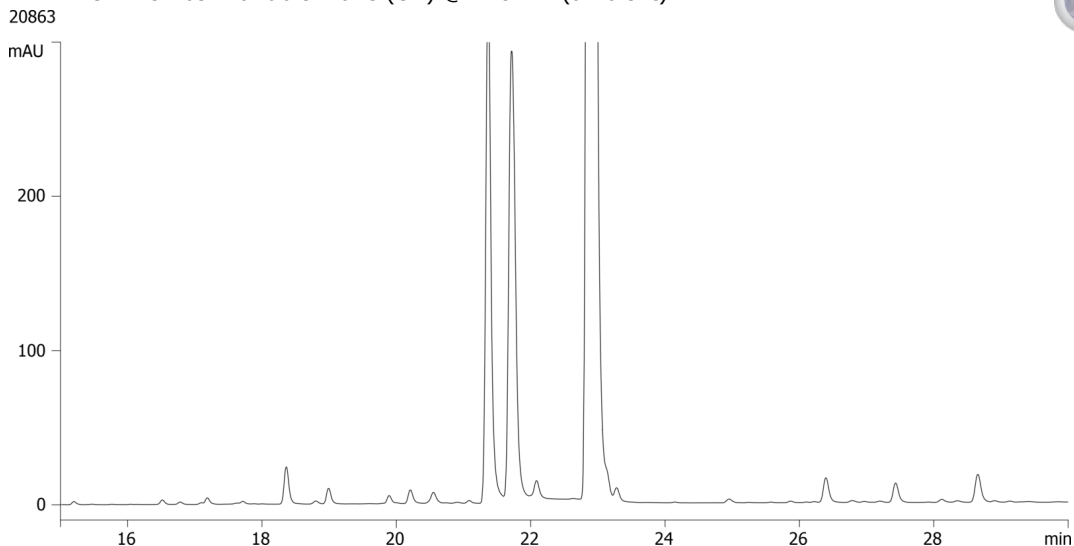