

Rapid analysis of BSA digest using Kinetex 1.3u C18 50x2.1mm

Column: Kinetex® 1.3 µm C18 100 Å, LC Column 50 x 2.1mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4515-AN

Elution Type: Gradient

Eluent A: Water with 0.1% TFA

Eluent B: Acetonitrile with 0.1% TFA


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1	95	5
	3	11	5	95
	4	11.5	5	95
	5	11.6	95	5

Flow Rate: 0.5 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ 210 nm (ambient)

21459

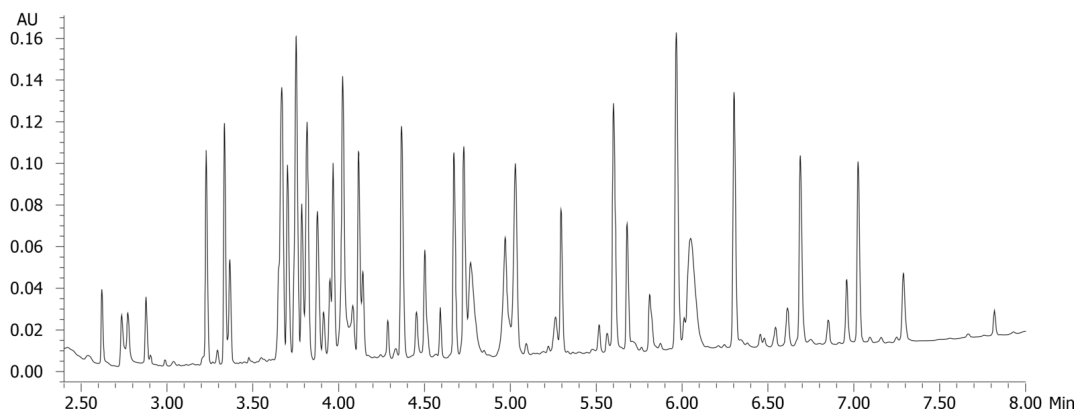

Products used in this application:



ANALYTES:

- 1 Bovine Serum Albumin (BSA)

