

Fast Pain Panel using Kinetex 2.6 u Biphenyl 50x3.0mm

Column: Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 3.0 mm, Ea

Dimensions: 50 x 3 mm ID

Order No: 00B-4622-Y0

Elution Type: Gradient

Eluent A: Water with 0.1% formic acid

Eluent B: Methanol with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	2.5	0	100
	3	3.5	0	100
	4	3.51	90	10
	5	5	90	10

Flow Rate: 0.7 mL/min


Col. Temp.: 40 °C

Detection: Tandem Mass Spec (MS-MS) @ (40 °C)

22121



Products used in this application:

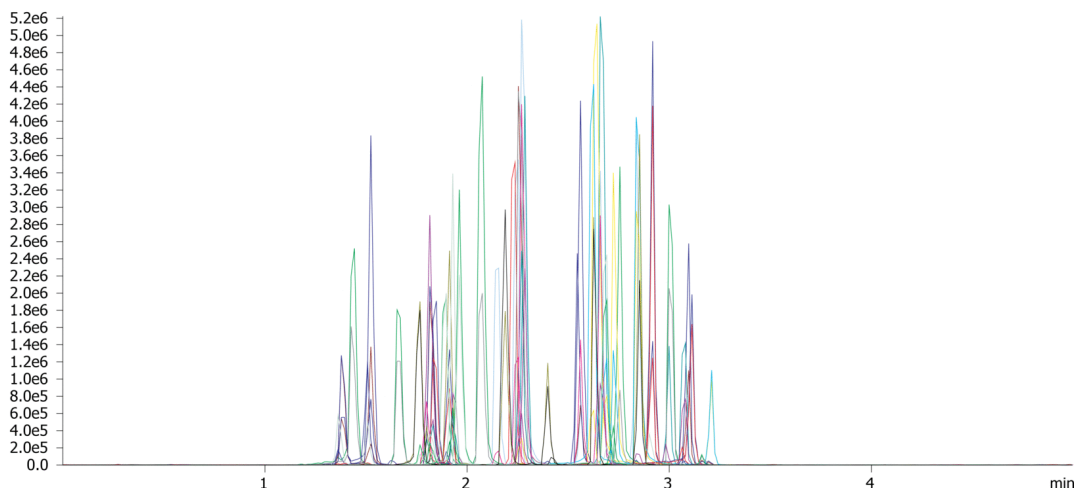

Fast Pain Panel using Kinetex 2.6 u Biphenyl 50x3.0mm

ANALYTES:

- | | | | |
|-----------|-------------------|-----------|----------|
| 1 | Morphine | 41 | Diazepam |
| 2 | Oxymorphone | | |
| 3 | Hydromorphone | | |
| 4 | Amphetamine | | |
| 5 | Naloxone | | |
| 6 | Methamphetamine | | |
| 7 | Codeine | | |
| 8 | MDA | | |
| 9 | Oxycodone | | |
| 10 | Naltrexone | | |
| 11 | Hydrocodone | | |
| 12 | MDMA | | |
| 13 | MDEA | | |
| 14 | Norfentanyl | | |
| 15 | Tramadol | | |
| 16 | Benzoyllecgonine | | |
| 17 | Meperidine | | |
| 18 | Meprobamate | | |
| 19 | Norbuprenorphine | | |
| 20 | Fentanyl | | |
| 21 | Buprenorphine | | |
| 22 | Flurazepam | | |
| 23 | Carisoprodol | | |
| 24 | PCP | | |
| 25 | Propoxyphene | | |
| 26 | Sufentanil | | |
| 27 | 6-MAM | | |
| 28 | Midazolam | | |
| 29 | Normeperidine | | |
| 30 | EDDP | | |
| 31 | Methadone | | |
| 32 | Lorazepam | | |
| 33 | Clonazepam | | |
| 34 | Norpropoxyphene | | |
| 35 | Oxazepam | | |
| 36 | Hydroxyalprazolam | | |
| 37 | Nordiazepam | | |
| 38 | Flunitrazepam | | |
| 39 | Temazepam | | |
| 40 | Alprazolam | | |

