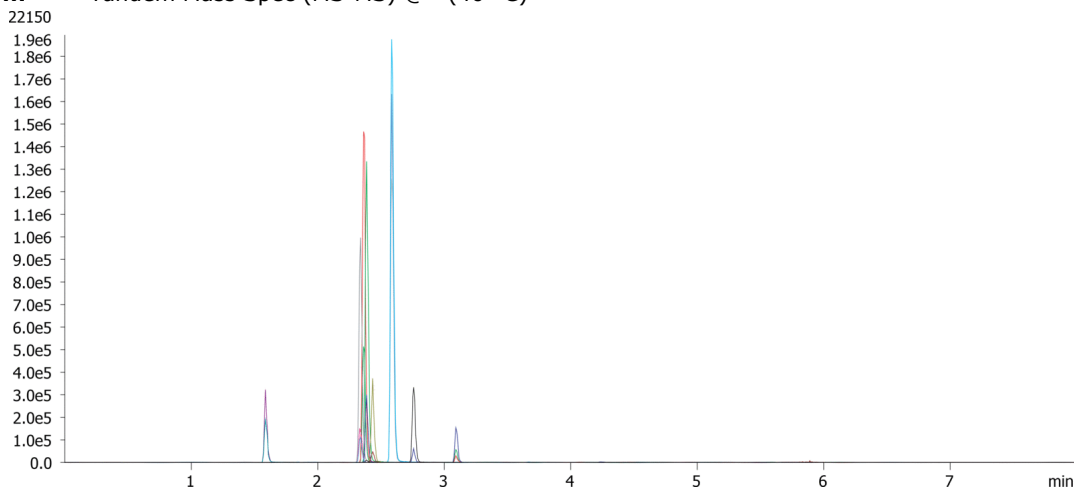


Antipsychotic Drugs by LCMS using Kinetex 2.6u Biphenyl 50x2.1mm**Column:** Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4622-AN**Elution Type:** Gradient**Eluent A:** 0.1% formic acid in water**Eluent B:** 0.1% formic acid in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	5	95
	3	5.1	95	5
	4	8	95	5

Flow Rate: 500 mL/min**Col. Temp.:** 40 °C**Detection:** Tandem Mass Spec (MS-MS) @ (40 °C)**ANALYTES:**

- 1** Olanzapine
Retention Time: 1.59 min
- 2** Clozapine
Retention Time: 2.34 min
- 3** 9-hydroxy Risperidone
Retention Time: 2.37 min
- 4** Risperidone
Retention Time: 2.39 min
- 5** Norquetiapine
Retention Time: 2.43 min
- 6** Quetiapine
Retention Time: 2.59 min
- 7** Ziprasidone
Retention Time: 2.76 min
- 8** Aripiprazole
Retention Time: 3.09 min

**Products used in this application:**