

## Batch Reproducibility with ABN 7 mix- Kinetex 2.6u F5 50x4.6mm

**Column:** Kinetex® F5 2.6µm 100 Å, LC Column 50 x 4.6 mm, Ea

**Dimensions:** 50 x 4.6 mm ID

**Order No:** 00B-4723-E0

**Elution Type:** Gradient

**Eluent A:** 0.1% Formic Acid in Water

**Eluent B:** 0.1% Formic Acid in Methanol


| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 95    | 5     |
|                   | 2        | 5          | 5     | 95    |

**Flow Rate:** 1.85 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

22803

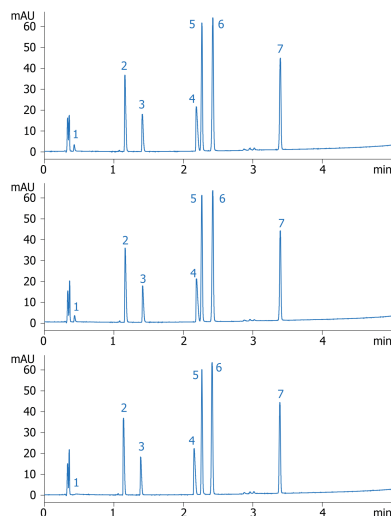

### ANALYTES:

- 1 Uracil (1mg/mL)
- 2 Pindolol (5mg/mL)
- 3 Chlorpheniramine (5mg/mL)
- 4 Nortriptyline (5mg/mL),
- 5 3-methy-4- nitrobenzoic acid (5mg/mL)
- 6 5-Methylsalicylaldehyde (5mg/mL)
- 7 Hexanophenone (20mg/mL)



Products used in this application:

