

## Stevia Standards on Kinetex 1.7u C18 150x2.1mm

**Column:** Kinetex® 1.7 µm C18 100 Å, LC Column 150 x 2.1 mm, Ea

**Dimensions:** 150 x 2.1 mm ID

**Order No:** 00F-4475-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% Formic Acid in Water

**Eluent B:** 0.1% Formic Acid in Acetonitrile


| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 72    | 28    |
|                   | 2        | 2          | 72    | 28    |
|                   | 3        | 2.01       | 68    | 32    |
|                   | 4        | 5          | 68    | 32    |
|                   | 5        | 12         | 1     | 99    |
|                   | 6        | 16         | 1     | 99    |
|                   | 7        | 17         | 72    | 28    |
|                   | 8        | 25         | 72    | 28    |

**Flow Rate:** 0.2 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 202 nm (ambient)

22913

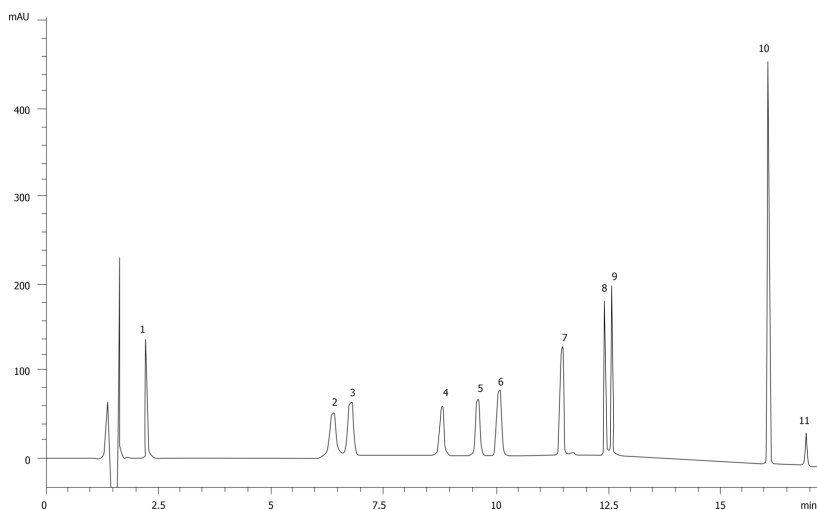

Products used in this application:



## Stevia Standards on Kinetex 1.7u C18 150x2.1mm

### ANALYTES:

- 1** Rebaudioside D
- 2** Rebaudioside A
- 3** Stevioside
- 4** Rebaudioside F
- 5** Rebaudioside C
- 6** Dulcoside A
- 7** Rubsoside
- 8** Rebaudioside B
- 9** Steviobioside
- 10** Steviol
- 11** Isosteviol

