

# HPLC Application

ID No.: 23056

## Fast Analysis of Sugars in Fruit Juices - Moscato (White Wine)

**Column:** Rezex<sup>™</sup> RCM-Monosaccharide Ca+2 (8%), LC Column 100 x 7.8 mm, Ea

**Dimensions:** 100 x 7.8 mm ID

**Order No:** 00D-0130-K0

**Elution Type:** Isocratic

**Eluent A:** Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

**Flow Rate:** 0.4 mL/min

**Col. Temp.:** 80 °C


**Detection:** Electrospray Mass Spec (ESMS) @ 0.000000000 (35 °C)  
23056

**REZEX**<sup>™</sup>  
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

**REZEX**<sup>™</sup>

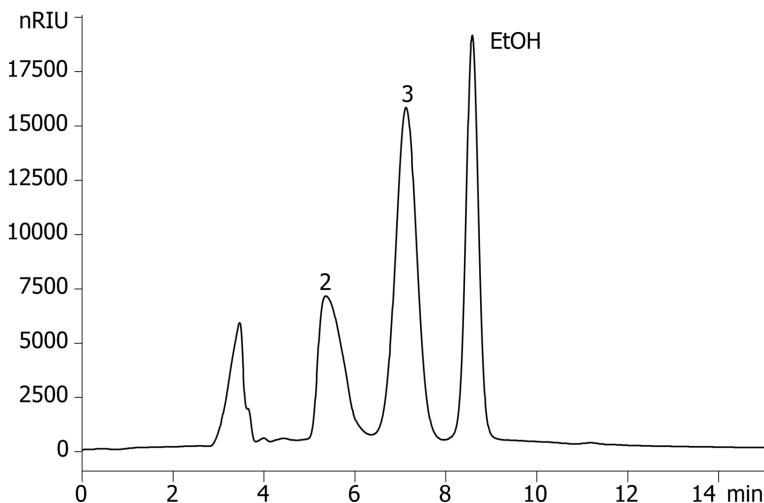

### ANALYTES:

- 1 Moscato (White Wine)

