

Fulvestrant and Fulvestrant enantiomer Lux Amylose 1 USP System Suitability


Column: Lux® 5 µm Amylose-1, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4732-E0
Elution Type: Isocratic
Eluent A: 2-methylpentane
Eluent B: Ethanol
Gradient Profile:

Step No.	Time (min)	Pct A	Pct B
1	0	88	12

Flow Rate: 1 mL/min
Col. Temp.: 40 °C
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (25 °C)



Products used in this application:

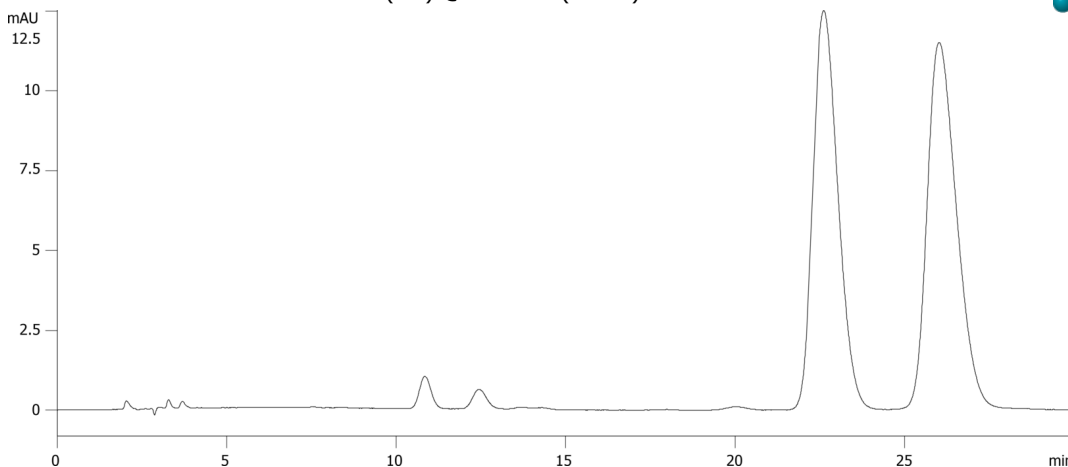

ANALYTES:

1 Fulvestrant enantiomer

Retention Time: 22.57 min

2 Fulvestrant

Retention Time: 25.578 min

