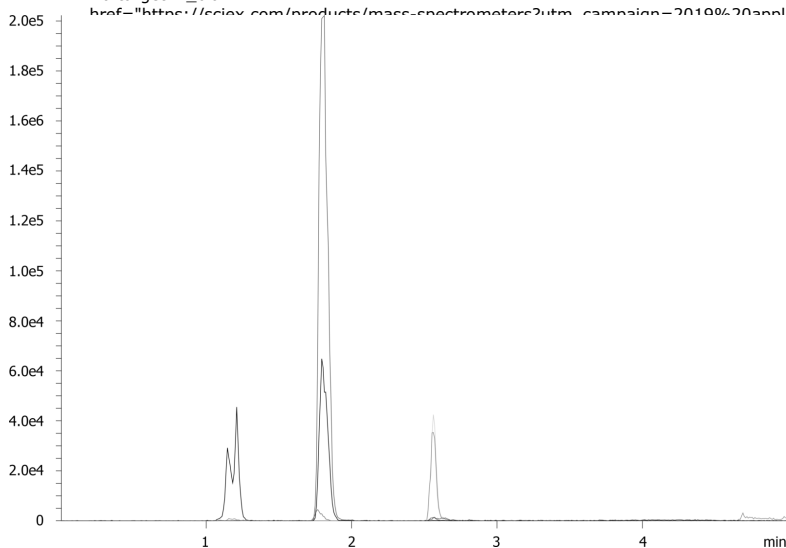


EtG/EtS in Urine using Synergi 2.5µm Hydro-RP, 100x3mm (3 min)**Column:** SynergiTM 2.5 µm Hydro-RP 100 Å, LC Column 100 x 3 mm, Ea**Dimensions:** 100 x 3 mm ID**Order No:** 00D-4387-Y0**Elution Type:** Gradient**Eluent A:** 10mM ammonium formate + 0.1% formic acid in water**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	1.5	5	95
	3	2	5	95
	4	2.01	100	0
	5	3	100	0

**Flow Rate:** 0.7 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (550 °C)**Detector Info:** SCIEX**ANALYTES:****1** Urine Interference

Retention Time: 0.73 min

2 Ethyl Sulfate (EtS)

Retention Time: 1.08 min

3 Ethyl Glucuronide (EtG)

Retention Time: 1.54 min

