


## Follitropin Reference Standard on Yarra SEC-X150, Impurity Profile

**Column:** N/A  
**Dimensions:** N/A  
**Order No:** N/A  
**Elution Type:** N/A  
**Eluent A:** 100 mM Sodium Phosphate, 200 mM Sodium Sulfate, pH 6.7  
**Gradient**  

| Step No. | Time (min) | Pct A |
|----------|------------|-------|
| 1        | 0          | 100   |

**Profile:**  
**Flow Rate:** 0.35 mL/min  
**Col. Temp.:** ambient  
**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 215 nm (ambient)  
 23793

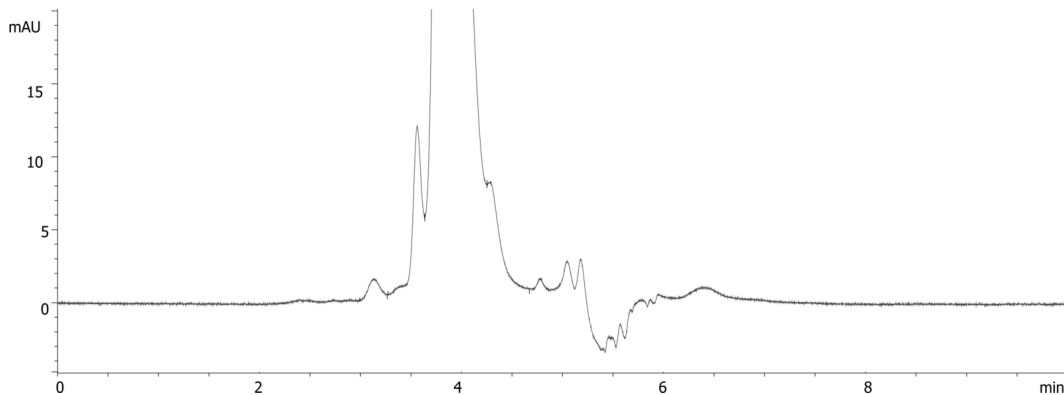

### ANALYTES:

1 Follitropin

