

Chiral Separation of 2-Chloromandelic Acid on Lux 3 μ m i-Cellulose-5**Column:** Lux® 3 μ m i-Cellulose-5, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4755-E0**Elution Type:** Isocratic**Eluent A:** 0.1% TFA in Dichloromethane**Eluent B:** 0.1% TFA in IPA**Gradient**


Step No.	Time (min)	Pct A	Pct B
1	0	99	1

Profile:**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 270 nm (ambient)

23880



Products used in this application:

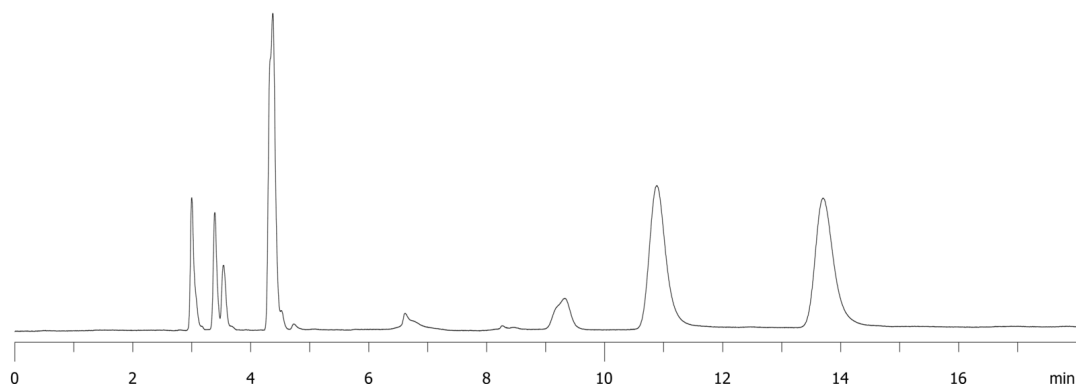

ANALYTES:

1 2-Chloromandelic acid

Retention Time: 10.881 min

2 2-Chloromandelic acid

Retention Time: 13.7 min

