

Chiral Separation of Sulpiride on Lux 3µm i-Cellulose-5


Column: Lux® 3 µm i-Cellulose-5, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4755-E0
Elution Type: Isocratic
Eluent A: 0.1% DEA in Hexane
Eluent B: 0.1% DEA in THF
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	60	40

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
 23888



Products used in this application:

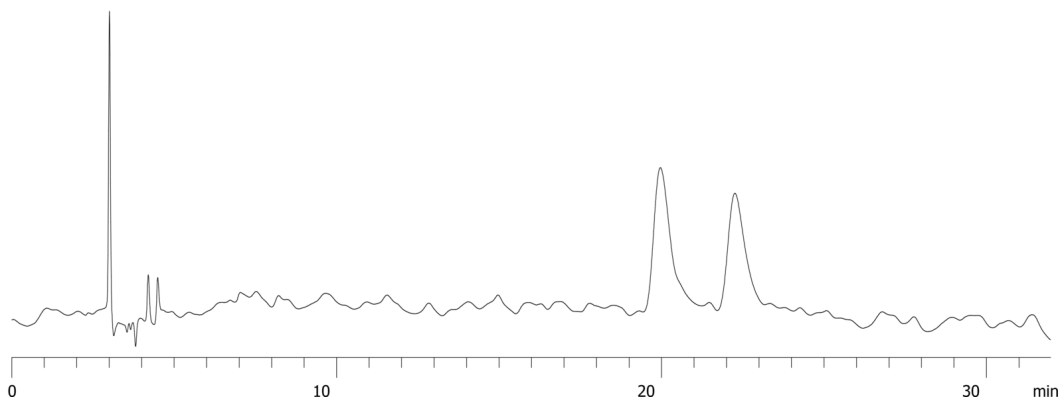

ANALYTES:

- 1 Sulpiride
Retention Time: 19.96 min
- 2 Sulpiride
Retention Time: 22.254 min

