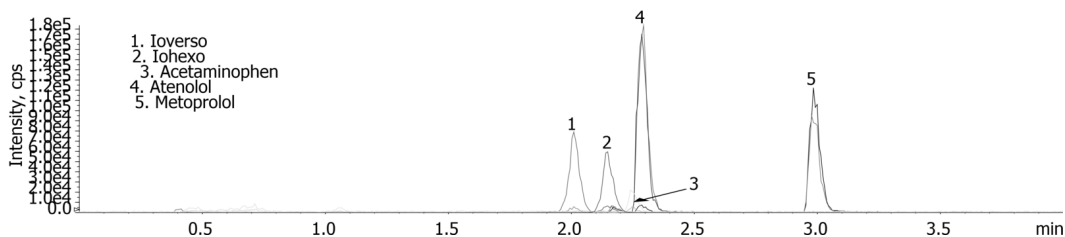


Iohexol in Serum - PPT (HCL04)**Column:** Kinetex® 2.6 µm Polar C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4759-AN**Elution Type:** Gradient**Eluent A:** 10 mM Ammonium Acetate (unadjusted pH)**Eluent B:** 10 mM Ammonium Acetate in MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0.01	100	0
	2	0.3	100	0
	3	2.5	35	65
	4	3	35	65
	5	3.1	5	95
	6	4.1	5	95
	7	4.5	100	0
	8	7	100	0

Flow Rate: 300 µL/min**Col. Temp.:** 40 °C**Detection:** Tandem Mass Spec (MS-MS) @ (600 °C)**Detector Info:** SCIEX

24250

**Products used in this application:****ANALYTES:****1** Ioversol

Retention Time: 2.01 min

2 Iohexol

Retention Time: 2.14 min

3 Acetaminophen

Retention Time: 2.24 min

4 Atenolol

Retention Time: 2.29 min

5 Metoprolol

Retention Time: 2.98 min

