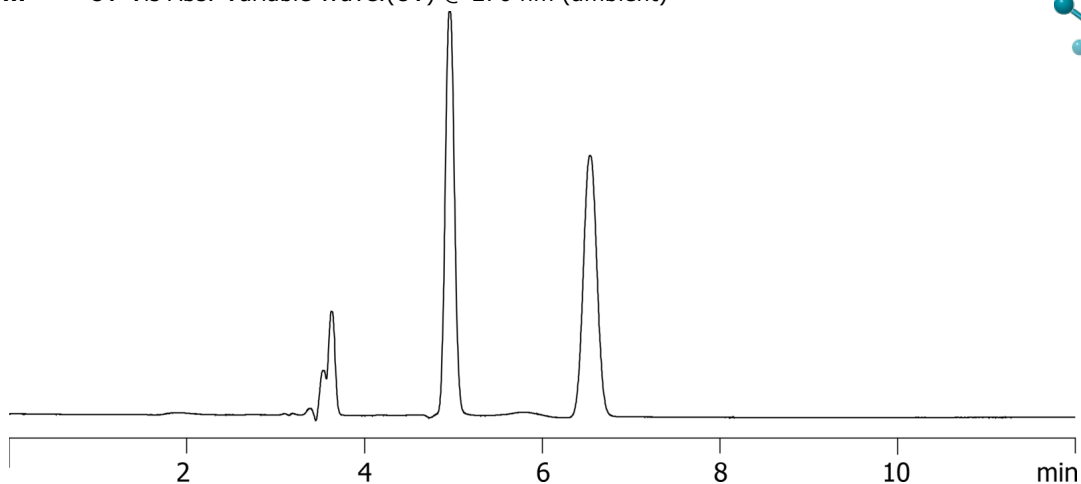


Chiral Separation of Tramadol on Lux 5 μ m i-Amylose-1 in NP**Column:** Lux[®] 5 μ m i-Amylose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4762-E0**Elution Type:** Isocratic**Eluent A:** 0.1% DEA in Hexane**Eluent B:** 0.1% DEA in IPA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	90	10

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 270 nm (ambient)**Products used in this application:****ANALYTES:****1** Tramadol

Retention Time: 4.958 min

2 Tramadol

Retention Time: 6.537 min

