


**Chiral Separation of Nadolol on Lux 5 $\mu$ m i-Amylose-1 in NP****Column:** Lux<sup>®</sup> 5  $\mu$ m i-Amylose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4762-E0**Elution Type:** Isocratic**Eluent A:** 0.1% DEA in Hexane**Eluent B:** 0.1% DEA in IPA

| Gradient | Step No. | Time (min) | Pct A | Pct B |
|----------|----------|------------|-------|-------|
| Profile: | 1        | 0          | 85    | 15    |

**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)**Products used in this application:****ANALYTES:**

- 1 Nadolol  
Retention Time: 12.355 min
- 2 Nadolol  
Retention Time: 12.915 min
- 3 Nadolol  
Retention Time: 22.312 min
- 4 Nadolol  
Retention Time: 26.359 min

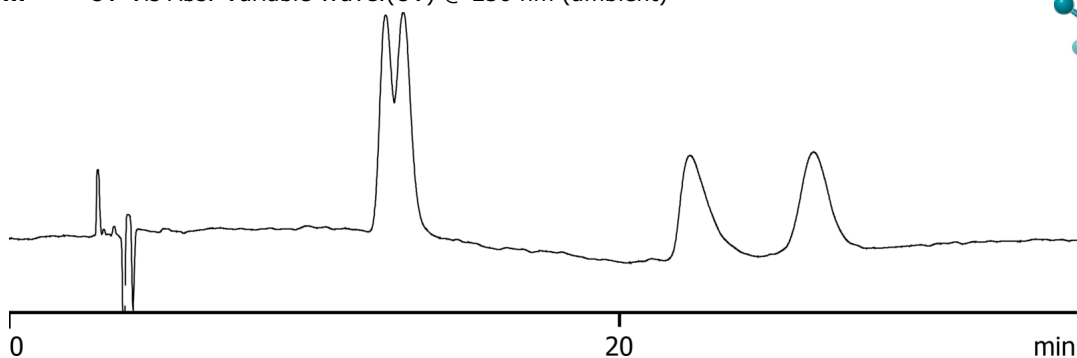