

Chiral Separation of Carbinoxamine on Lux 5µm i-Amylose-1 in NP


Column: Lux® 5 µm i-Amylose-1, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4762-E0
Elution Type: Isocratic
Eluent A: 0.1% DEA in Hexane
Eluent B: 0.1% DEA in THF
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	90	10

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: Electrospray Mass Spec (ESMS) @ 0.000000000 (ambient)
24505



Products used in this application:

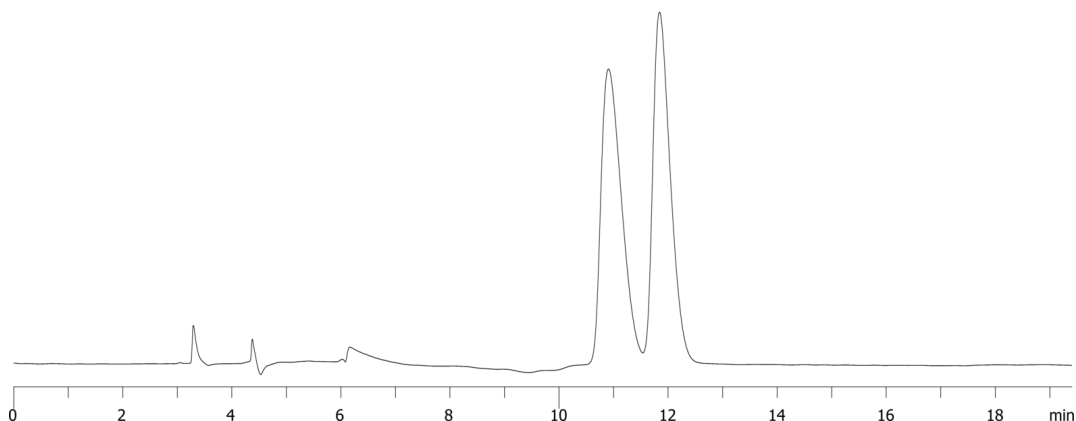

ANALYTES:

- 1 Carbinoxamine
Retention Time: 10.904 min
- 2 Carbinoxamine
Retention Time: 11.84 min